


80K Body Sculpting Machine

< Instruction Manual >



To ensure proper use of this product, please read the instruction manual carefully and store it in a safe place for future reference.

Accessories Details

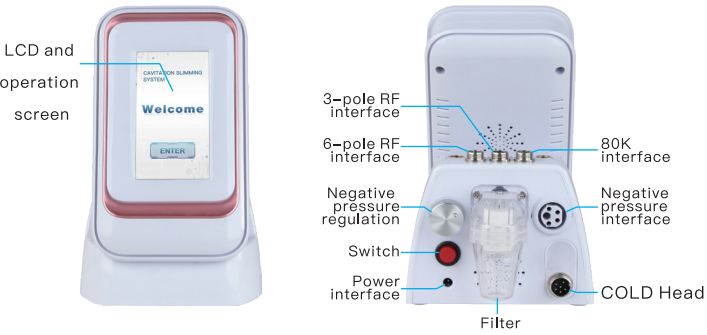


Function of Different Head

3 RF Head: Suggest Use Energy Level 1–2

For Face and Eye Use


3 RF head firms and lifts the skin, reduces wrinkles, and improves skin texture. It also enhances eye circulation, shrinks pores, and reduces puffiness. Additionally, it helps eliminate wrinkles, dark circles, and under-eye bags, while rejuvenating the skin and improving the contour of the eye area.

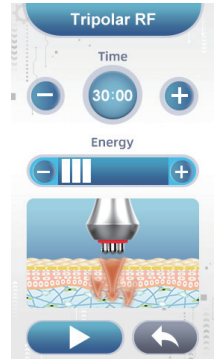

Ice pack freezer head

For Face and Neck Use


After using radiofrequency on the face, it can effectively prevent the continued rupture of capillaries, and ice compress can cause local blood vessel contraction, reduce blood circulation, and inhibit inflammatory reactions.

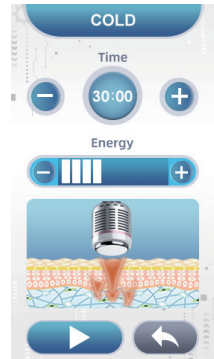

6 RF Head: Suggest Use Energy Level 2–5

For Back, Waist, Thighs, and Calves Use

6 RF head helps tighten the waist, firm the body, and enhance overall body shaping. It also softens the muscles of the arms and thighs, provides a full-body massage, and tones the arms and thighs for a more sculpted appearance.

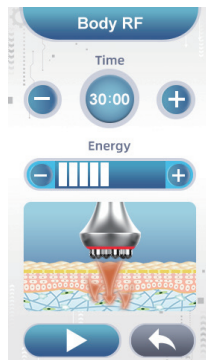

Negative Pressure Massage Head

Body Maintenance

With diversified level adjustment, negative pressure massage ensures a comfortable and personalized body care experience. It improves skin texture by stimulating collagen regeneration and enhancing elasticity, while also addressing skin sagging caused by fat loss and aging.




Additionally, it helps eliminate waste materials from the body and boosts metabolism. The negative pressure massage function promotes fat cell reduction, while the vacuum negative pressure massage effectively absorbs various skin tissues. Combined with multi-polar thermal melting radiofrequency, it allows the radiofrequency energy to penetrate deeper into the skin tissue more efficiently. At the same time, the negative pressure adsorption enhances blood circulation and metabolism, fully relaxing muscle tissue and relieving skin and muscle fatigue.

80K Fat Burn Head

For Body Shaping Use

Fat burn head helps dissolve fat, detoxify the lymphatic system, tighten the skin, and enhance skin elasticity. It directly targets the fat layer with directional RF (radiofrequency) output, accelerating fat cell movement to generate heat energy. This process promotes the elimination of excess fat and toxins from the body through sweat glands, liver–intestine circulation, and lymphatic circulation, effectively aiding in fat reduction.

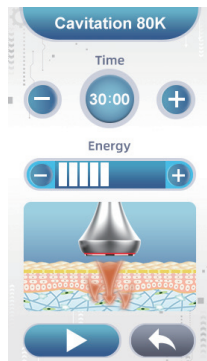

Operation Instructions

1. Connect the hand piece to the corresponding jack and connect the power cord.
2. Press the power switch, the indicator light and LCD will be on.
3. Choose the function you need.
4. Touch the 80K button to display figure 2, adjust the required mode, working time, working intensity, press the “▶” button, the 80K burn head starts to work, and other R.F head function selection methods are the same operation.
5. After using, press the “◀” button to return the interface shown in Figure 1, and then turn off the power button.



P1




P2

Using Guidance


Thigh Treatment

- (1) Push upward from the bottom of the thigh toward the groin to promote lymphatic drainage.
- (2) Use circular motions in a counterclockwise direction, moving from the bottom to the top, to help break down fat.
- (3) Push back and forth to target and decompose fatty tissue.
- (4) Gently pull downward from the knee to the groin to enhance and refine the thigh curve.

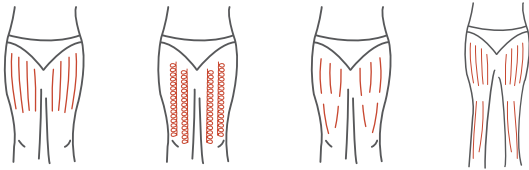


Hips Treatment

- (1) Gently pull upward toward the waist, following the direction of lymphatic flow.
- (2) Use circular motions in a counterclockwise direction, moving from the bottom to the top, and then glide back and forth.
- (3) Following the natural muscle texture, glide and compact the skin to enhance the buttocks' fullness and firmness, creating a perfect, sculpted curve.

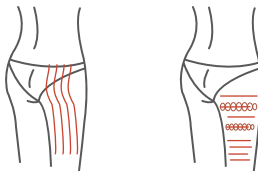

Back Treatment

- (1)Glide back and forth 2–3 times to relax the muscles.
- (2)Use circular motions in a counterclockwise direction to stimulate the lymph nodes.
- (3)Move in a counterclockwise motion while gliding up and down the back to help break down fat.
- (4)Apply lymphatic drainage techniques to guide toxins toward the lymph nodes, improving the back's contour and curve.

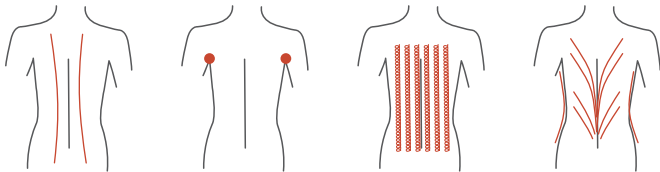

Arm Treatment

- (1)Glide from the inner elbow toward the armpit to perform lymphatic drainage.
- (2)Use circular motions to stimulate the lymph nodes.
- (3)Break down excess fat by applying counterclockwise circular motions or gliding back and forth.
- (4)Firm and tighten the skin by gliding upward from the elbow to the armpit.



Face Treatment

(1) Apply the product evenly to the face. Using the RF head, glide from the bottom of the face upward to lift and tighten the skin. Maintain a gentle tension on the skin during the process. It is recommended to use for 5 minutes.


(2) Apply the product evenly to the face. Using the RF head, follow the direction of the arrows as indicated show. It is recommended to use for 5–8 minutes.



Abdomen Treatment

(1) Starting from the navel, use clockwise circular motions, gradually expanding from small to large circles. This helps stimulate peristalsis in the large intestine.

(2) Transition from small to large circles while pushing in a counterclockwise direction to further enhance the treatment effect.

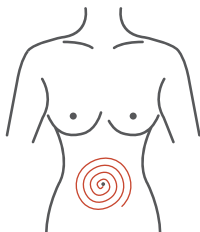

Warm Tips for You

1. Do not use on dry or bare skin, as this may cause skin burns. Always apply cream or oil to the skin before use.
2. Ensure an aqueous medium (e.g., essential oil or cream) is applied during use. Operating the device without a medium may damage the machine. Avoid using on bony areas and do not linger on the same body part for too long.
3. Do not connect the power until the probe is properly installed. Always use an aqueous medium (e.g., essential oil or cream) to prevent damage to the machine.
4. Do not replace the probe while the machine is charging.
5. Remove all metal jewelry before using the machine to avoid interference or injury.
6. Allow the machine to rest after continuous use for 1 hour. Stop using it for at least 10 minutes to prevent overheating.
7. Turn off the power after each use and clean the probe thoroughly. Proper maintenance ensures better performance for future use.
8. Keep the probe clean after every use. Wipe it with a dry cloth—do not use water to clean the probe.
9. After treatment, cleanse your skin and apply a hydrating mask, cream, or oil after 10 minutes to soothe and nourish the skin.
10. When using the 80K Head, some users may experience tinnitus. This is normal, as the sound waves produced may exceed the range of human hearing.

Warranty & After-Sales Service

We guarantee that every machine undergoes rigorous quality checks to ensure it meets the required standards. We offer a 1-year warranty, and if you encounter any technical difficulties during use, our professional team will provide expert assistance to promptly resolve your issues.

We sincerely hope that you will join us on this journey and enjoy the endless experiences our products have to offer.

If you have any questions or concerns, please do not hesitate to contact us. We are always here to support you.

