





USER MANUAL

► INTRODUCTION

Thank you for purchasing the LED red light therapy hat, which may be the most advanced phototherapy solution in the market. We have not only made them effective, but also had many user testings in order to provide the best experience. We are very happy that you have the new generation of IR light therapy hats, which has everything you need for safety, effectiveness, portability, and care. These light wavelengths have been scientifically proven to be helpful when used. Please read this manual thoroughly before you begin using this product in order to get the best effects.


NOTICE

The product cannot replace medicine, they can only serve as auxiliary functions.

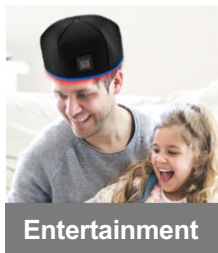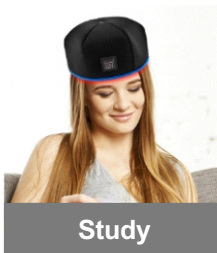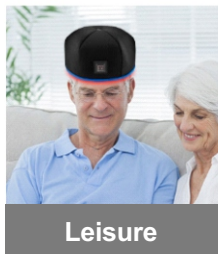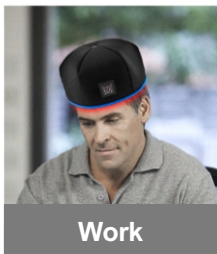


Multiple use cases

▶ PRODUCT DIAGRAM

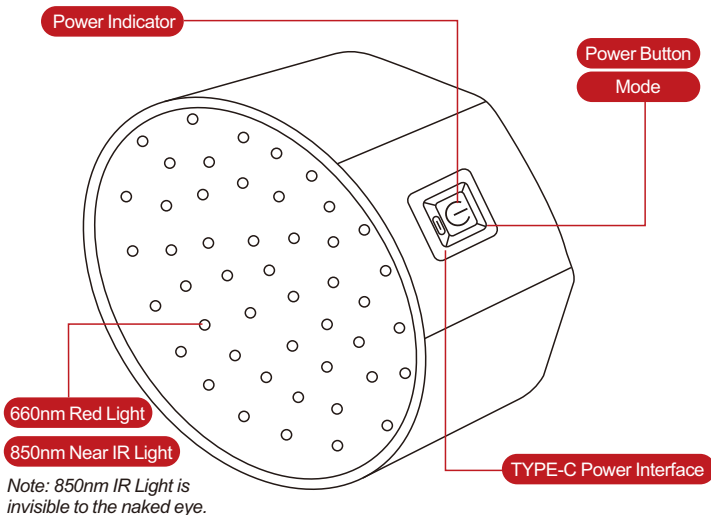

WHAT'S INCLUDED

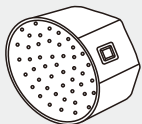

- Cap x 1



- Use manual x 1



- TYPE-C cable x 1



- Packaging box x 1

▶ TWO STYLES OF LED RED LIGHT THERAPY CAP

Three modes with wired

Normal Mode

Red light and IR light work simultaneously, penetrating deep into the scalp, stimulating blood circulation in hair follicles, improving cell metabolism, helping hair health, increasing hair density and health.



The indicator light shows red

10Hz Pulse Mode

Flashing 10 times per second increases dopamine levels in the brain, improves mood and cognitive abilities.



The indicator light shows green

40Hz Pulse Mode

Flashing 40 times per second, stimulating the cerebral cortex, accelerating metabolism, promoting photobiological regulation, and helping hair health.



The indicator light shows blue

Three modes with wireless

Normal Mode

Red light and IR light work simultaneously, penetrating deep into the scalp, stimulating blood circulation in hair follicles, improving cell metabolism, helping hair health, increasing hair density and health.



Long press to blue light

10Hz Pulse Mode

Flashing 10 times per second increases dopamine levels in the brain, improves mood and cognitive abilities.



No changes in blue light


40Hz Pulse Mode

Flashing 40 times per second, stimulating the cerebral cortex, accelerating metabolism, promoting photobiological regulation, and helping hair health.



No changes in blue light

1. Long press to turn off
2. Charging - flashing red
3. Charging complete - solid red
4. Automatically shut off after 20 minutes

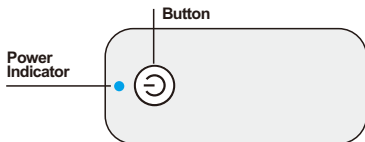


Controller

► HOW TO USE


1

Insert power adapter to wall outlet and connect the power cord to power switch on the cap.

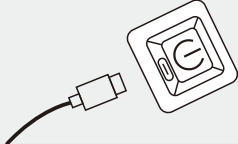


2

Long press the button to start the device, and lightly press the button to cycle through mode settings: light remains, pulsed light (10Hz), and pulsed light (40Hz). The naked eyes can only see red light, while IR light is invisible to the naked eyes.

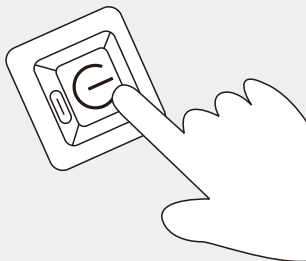

3

Put the Red Light Therapy Cap on your head.

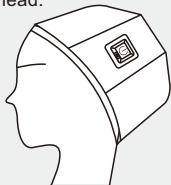

4


The usage time is 20 minutes, and it will automatically shut off after 20 minutes.



5

Long press 2 seconds to turn on / off.

2
seconds

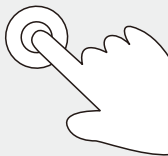

► THE FEATURES AND BENEFITS

Hair health is a long process that requires your attention and investment. It requires time, patience, and perseverance with daily use. Hundreds of clinical studies have proven that phototherapy has a positive impact on health, and we believe your body is benefiting from every use. But there are many factors that can affect the speed at which your symptoms improve.

Overall, the progression of health or skin related issues may require 8-12 weeks of consistent use. By consistently using the product every day, you can maximize the activity of your body cells.



If you have any questions about photobiological regulation, we suggest that you consult a trusted medical care provider or doctor.

Product Specifications

Product Name : LED Red Light Therapy Cap

Led Qty : 120pcs (dual chips, 120pcs 660nm chips+120pcs 850nm chips)

Wavelength : 660nm/850nm

Timer : 20 minutes

Material : Lycra

Diameter : 20 cm

Weight : 138 g

Power : Type-c 5V-2A

Battery Capacity : 1800MHA (wireless)

► WARNING

- 1 The red light+infrared therapy lamp is an electronic device.
To avoid electric shock : * Do not stand on the device.**

- * Do not wash electric parts.
- * Do not touch plugs or devices with wet hands.
- * Do not use the device when cords or plugs are damaged.
- * Do not use the device if it is damaged or comes in contact with water.
- * Do not carry or lift up the device by the power cord.
- * Do not use the power cord as a handle to prevent it from being caught by the door, or pull the power cord towards sharp edges or corners.
- * Do not use the device near hot surfaces or environment.
- * Do not fold and place heavy objects on the device.
- * Store the device with the original packaging to avoid any kind of damage.

- 2 Keep the product and its components out of reach of children.**

- 3 Do not use this device on open wounds as high temperatures may increase circulation and cause bleeding.**

- 4 Do not use with liniments, ointments, or ointments containing febrile ingredients, as it may cause skin burns.**

- 5 Unsuitable for pregnant and photosensitive individuals**

To check if you have photosensitive skin, hold your arms 6 inches away from the device for 3 minutes. If redness or discomfort occurs, you may have photosensitive skin. In this case, please stop using and consult a doctor.

► WARRANTY

This product provides a one-year free warranty and lifetime maintenance service. Within one year from the date of purchase, if there is a product problem (non-artificial damage), please contact us and we will be happy to serve you. Note : This Warranty does not cover any damage to this Product which results from accident, abuse, misuse or non-compliance with the user manual.

THANKS