



Facial + Neck Mask

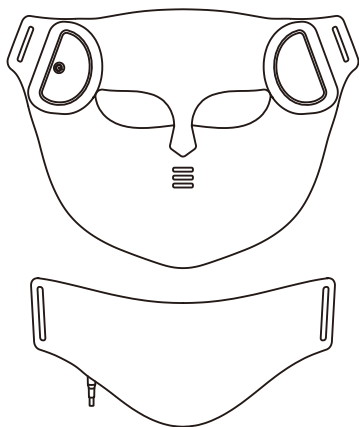

- operation instructions -

PRODUCT INTRODUCTION

The beauty mask utilizes advanced focused ultra-narrow light technology, boasting an energy efficiency six times higher than traditional LED light sources. Through four specific wavelengths of light, 100 LED beads, and 400 light sources, it precisely targets different layers of the skin, achieving comprehensive skin care. The mask is designed to fit the contours of the face, making it convenient to use, allowing you to enjoy a professional-level beauty and skincare experience at home.

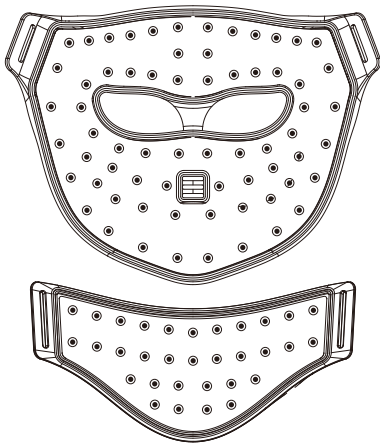






Suggest using 2-3 times a week for better results

PREFACE

Thank you for purchasing our product! With this LED face & neck mask instrument, you or your loved ones can enjoy a skincare experience for younger, clearer and brighter skin.

For a better user experience and results, please read this manual carefully before use. Also, upon receiving the product, check for damage and ensure all accessories are complete.

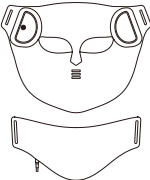INVENTORY

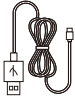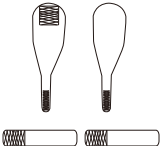
<p>1</p> 	<p>2</p> 	<p>3</p> 
<p>4</p> 	<p>5</p> 	
<p>① Beauty instrument x2 ② Instruction manual x1 ③ Controller x1</p> <p>④ Charging cable x1 ⑤ Bandage x 4</p>		

FIRST TIME SETUP

Before opening the package to use, please be sure to give the mask controller handle fully charged. When charging is in progress, the battery level display will flash in increments. Three lights on simultaneously indicate charging completion. Charging time is 2-3 hours.

The face & neck mask is equipped with adjustable straps, allowing you to wear it comfortably. When installing the face & neck mask, thread the strap through the entrance hole on the other side of the mask, wrap the strap around the head and secure it firmly.

WEARABLE CONNECTION


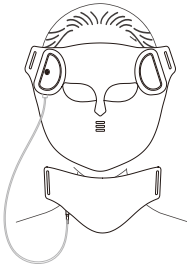


After the photonics device is installed, adjust the strap buckle to a comfortable state, Connect to the controller and turn on the photon.


Before use, please illuminate the face for 1 minute, and continue to use it without discomfort symptoms.

Start from gear 1-2, establish skin tolerance and select energy gear according to skin actual conditions. Do not blindly increase the time and frequency of light, and adjust the light time according to the actual situation of the skin. It is recommended that this product be used 3-4 times a week, and a single use time is 15-25 minutes


CHARGING METHOD

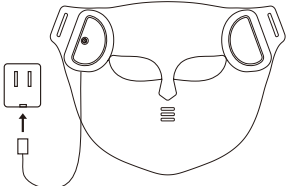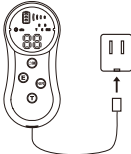Host charging method	Remote control charging method
 <p>⚠ Please fully charge the device when using it for the first time. Charging indicator light red, fully charged indicator light green</p>	 <p>When charging is in progress, the battery level display will flash in increments, (From one grid to three grids, a total of three grids), filled with three grids that are constantly bright.</p>

PHOTON INTERPRETATION

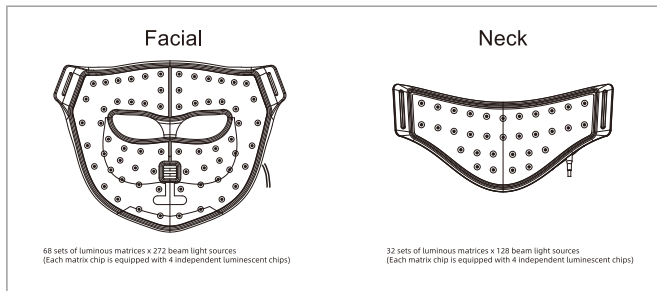

INSTALLATION DIAGRAM

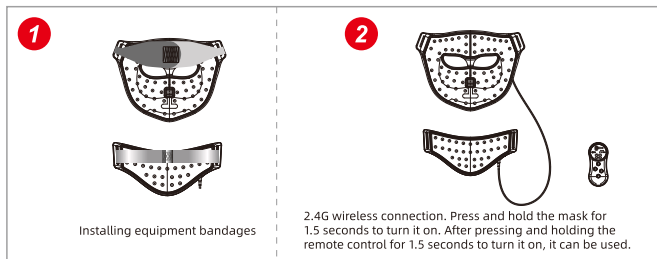

Tips

Before use, please irradiate the face for 1 minute, no discomfort symptoms can continue to use.

Please start using from low gear, build skin tolerance and choose the energy level according to the actual skin condition. Please do not blindly increase the time and frequency of light, you need to adjust the length of light according to the actual skin condition. It is recommended to use this product 3 - 5 times a week, a single use time of 10 minutes, a continuous cycle of 4 weeks.

CONTROLLER OPERATION EXPLANATION


Host



Long press to turn on/off (it will automatically turn off after 10 minutes of operation by default),
Short press to switch light sources, with seven types of light sources available for cyclic switching.

Power on/off/
Mode switching

Charging/
plugging port

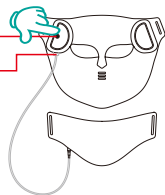

Remote control



!Please turn on the host switch first! When the host is turned on, long press the power button to turn it on,
Short press to switch light sources, seven types of light sources can be cycled.




After the device is turned on, the infrared mode is turned off by default. Lightly press to switch infrared mode on/off.
When this mode is enabled,
The "NIR" light on the display screen is on.



After turning on the device, the default energy level is level 1,
Lightly press once to switch between enhanced and 3-speed energy cycle adjustment.




After the device is turned on, the default time is 10 minutes, which is displayed as "10" on the screen.
Lightly press the T key to increase the time by 5 minutes.
10/15/20/25/30 minutes, 5 levels can be cyclically adjusted.

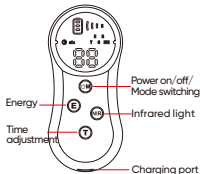


- ✓ When the device is in Red Light mode, short press M to switch to Blue Light mode. The screen will display "B", all blue lights on the FPC panel will illuminate, the timer will show the remaining countdown time, and the energy level will be the default/customized setting.
- ✓ When the device is in Blue Light mode, short press M to switch to Green Light mode. The screen will display "G", all green lights on the FPC panel will illuminate, the timer will show the remaining countdown time, and the energy level will be the default/customized setting.
- ✓ When the device is in Green Light mode, short press M to switch to Yellow Light mode. The screen will display "Y", all yellow lights on the FPC panel will illuminate, the timer will show the remaining countdown time, and the energy level will be the default/customized setting.
- ✓ When the device is in Yellow Light mode, short press M to switch to Purple Light mode. The screen will display "RB", both red and blue lights on the FPC panel will illuminate, the timer will show the remaining countdown time, and the energy level will be the default/customized setting.
- ✓ When the device is in Purple Light mode, short press M to switch to Cyan Light mode. The screen will display "BG", both blue and green lights on the FPC panel will illuminate, the timer will show the remaining countdown time, and the energy level will be the default/customized setting.
- ✓ When the device is in Cyan Light mode, short press M to switch to White light mode. The screen will display "RGB", red, blue, and green lights on the FPC panel will all illuminate. The timer will show the remaining countdown time, and the energy level will be the default/customized setting. Pressing M again will cycle through the modes sequentially.

SPECTRAL ANALYSIS

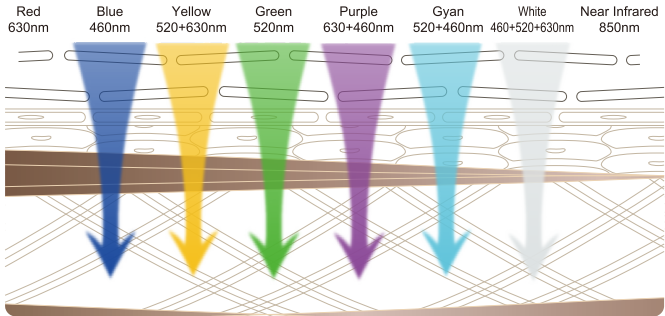

A study by NASA in the 20th century found that different wavelengths of light can have different skincare effects on the skin.

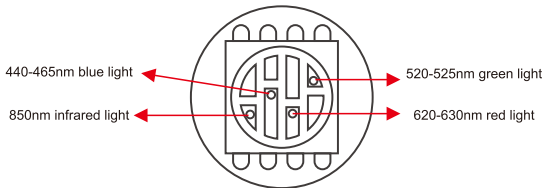

8 types of light waves effectively improve various skin problems:

- Red light (620 - 630nm) : Promotes collagen production, improves skin vitality, and reduces wrinkles and fine lines.
- Blue light (440 - 465nm) : Suitable for sensitive skin; reduces acne breakouts and improves acne - prone skin.
- Green light (520 - 525nm) : Improves sebum balance and regulates skin's water - oil equilibrium.
- Yellow light (630nm + 520nm) : Brightens skin tone, whitens skin, decomposes pigments, and improves scars.
- Purple light (630nm + 460nm) : Dual - wavelength light of red and blue; improves acne scars.
- Cyan light (525nm + 460nm) : Boosts skin energy and promotes metabolism.
- White light (460nm + 520nm + 630nm) : Penetrates deep skin layers for improvement.
- Near - infrared light (850nm) : Repairs skin and accelerates metabolic processes.

SAFETY GUIDELINES

Before using this product for the first time, please conduct a photosensitivity test

1. Place the device on your arm, turn on the switch, and let it emit light for 10 minutes. If symptoms such as skin redness or itching occur, you may have photosensitive substances. We suggest that you stop using the product or consult a beautician before using it again.
2. Before use, please clean your face and keep it dry (without makeup). Cosmetics may contain photosensitive ingredients that may affect the skin care effect.
3. After use, the skin may appear red, which is a normal reaction of the skin and will return to normal within 1 - 2 hours. If there is a burning sensation or any other discomfort during use, please stop using immediately.
4. Please apply a facial mask to replenish water after use and avoid direct sunlight. Repeat 3 times a week for 4 weeks for better results.

PRODUCT PARAMETERS:

Product name	Facial + Neck Mask		
Time timing	10/15/20/25/30 minutes	Net weight	365g
Charging input	5V/0.8A	Output power	3W
Main materials	silicone	Battery capacity	2000 mA
Charging duration	2-3 H		

