


英文一份

折页双面印刷

尺寸：80*140mm

材质：128克铜版纸

LED micro power beauty mask

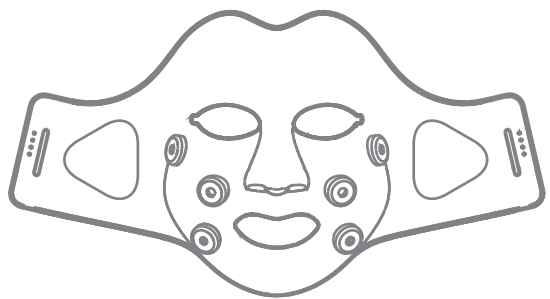

Read all the instruction before use

Product description:

LED micro power beauty mask is an Cross-age beauty product. This product uses micro-electronics in the face to massage acupuncture points and combines with LED specific beauty red light, 15 minutes one time, Reach the effect of facial shaping and beauty. When using this product, it can be used with facial masks or beauty products to promote absorption, improve blood circulation and brighten the complexion.

How to use:

1. First take out the beauty mask and install the micro power massage head on the beauty mask.

2. Wash and dry the face, apply a mask or skin care product to the face, and then adjust and fix the beauty mask on the face with straps.

3. Choose the mode according to your preferences, each mode is recommended for 5 minutes, (import with vibration mode. import with red led light and vibration mode. red led light beauty mode.)

4. After use, remove the massage head, take off the beauty mask, rinse with water and dry the mask.

5. First use, please control within 5 minutes, until the face to adapt to vibration massage, and then extend the use of time.

6. The face mask is recommended for using 2-3 times a week.

LED micro-power face shaping instrument operating instructions:

• Press “0n” key to turn on , “0FF” key to turn off ;

• MODE key has three modes, 1. Import with vibration mode;

2. Import with red led light and vibration mode; 3.Red led light beauty mode, The three modes are switched cycling, after turning on there is defaulted to vibration mode.

☉ Press this key to control vibration intensity, a total of three,

followed by high-frequency, intermediate-frequency, low-frequency, cycle switching.

☉ Press this button to select the motor for normal or cycle rotation; cycle rotation, only two of the three groups of motor at the same time, the interval of about 10 seconds automatically switch to the other two groups, its combination form is 1-2,2-3, 3-1. In this mode, only two groups are rotating.

• In vibration mode, it will automatically shut down after 5 minutes without any operation , LED lamp mode automatically shuts down after 15 minutes without any operation. Vibration plus LED light mode automatically shuts down after 5 minutes without any operation. Shutdown means that the microcontroller enters sleep mode.

Product Benefits:

1. Massage thin face effect: effective acupoint vibration massage for the face, micro-vibration frequency per minute 8000-12000 times, adhere to use of tight skin, v face shape, make the skin lightening.

2. Led Nano Red light effect: red light wavelength 630μm plus or minus 15μm, in skin care, health care effect is significant,known as active light. Red Light can improve the activity of cells, promote cell metabolism, so that the skin a large number of Collagen secretion and fibrous tissue to fill itself. Accelerate blood circulation, increase skin elasticity, improve the condition of dull skin, thus achieve anti-oxidation, repair the function of damaged cells.

3.Germanium titanium energy element effect: our company produced micro-power cosmetic mask added Germanium titanium and other energy elements, using the semiconductor function of energy metal, promote blood circulation, and through the contact with human skin up anti-fatigue health care effect.

4.Mask and skin care products auxiliary effect:

The traditional mask in about 5 minutes began to dry up, 40% of the essence is evaporated by the air instead of being absorbed by the face . Face mask or skin care products combined with micro-powered beauty mask using, not only can prevent the loss of mask by air evaporation, but also can make mask and skin care products appear to increase the effect.

parameter:

Product name: LED micro power beauty mask

Model: TM-MR-097

Weight: 227g

Voltage: 4.5V

Power: <0.975W

Red light wavelength:630μm plus or minus 15μm

Specification: 195(L)*400(w)*0.2(thickness) mm

Power supply: 3 AA alkaline batteries

Material: silicone + electronic components + ABS

Notes:

1. Please read the notes carefully before using this product.

2. Please use the weak vibration intensity first when the first time, the time does not exceed 5 minutes.

3. Do not use for purposes other than facial massage.

4. Do not use this product if there is a wound on the face.

5. When the red light flashes, it means the battery needs to be replaced. Please replace all batteries at one time. If not used for a long time, please take out the battery.

Warranty card

1. This warranty card free warranty is limited to the product due to production, technical and quality problems caused by the normal use of the product;

2. The warranty period shall be-year, that is, the date of sale of the product shall take effect, and shall be valid during the year;

3. During the warranty period, the company will replace the faulty parts free of charge.

4. If artificially damaged, the customer will have to pay the cost of maintenance spare parts and round trip freight.

5. The warranty card should be filled in block letters, the handwriting is clear and must not be altered. Please keep it properly. If you lose it, you will not make up for it.

product name		Product number	
user name		contact number	
Purchase date		Ways to purchase	
User address			
Maintenance records			
date	Maintenance content	Maintenance results	principal