Multi-functional Ultrasonic Beauty Device

(Rechargeable)

User Manual

Thank you for purchasing our product. Please ensure you read this manual thoroughly before using the device.


Precautions

1. This device is not suitable for individuals (including children) with diminished sensory or mental capacities, or those lacking experience and knowledge.
2. After using the device, please clean it promptly to prevent the buildup of oils or gel residues.
3. The internal battery of this device is not user-replaceable.
4. Exercise caution regarding water exposure when using this device to avoid risks of electric shock and other potential hazards.
5. This device is not washable. After use, wipe it clean with a damp cloth and then dry it with a dry towel.
6. Should you have any questions during the use of this product, please consult the {Frequently Asked Questions} section in the manual.

Prohibitions

1. Do not disassemble or modify the device on your own as it may cause malfunctions or pose dangers.
2. The device should not be used by children, pregnant women, or individuals with restricted physical mobility.
3. Avoid using damaged chargers or plugs that do not fit securely, as this can lead to electric shocks, short circuits, or fires.
4. If the device overheats during use or charging, disconnect the power immediately.
5. Avoid dropping or subjecting the device to impacts, which could lead to malfunctions.
6. Do not use this device if you are experiencing fever or any other physical abnormalities.
7. Avoid using this device if you have any of the following conditions: heart disease, hemophilia, allergic predispositions, or dermatitis.
8. Do not use the device on affected areas such as wounds, post-surgical sites, inside the mouth, upper eyelids, eyeballs, areas around the ears, or near the thyroid gland.

Product Structure

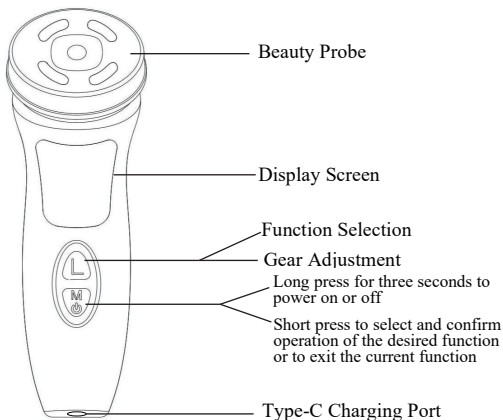Packing List

Ultrasonic Beauty Device *1

User Manual *1

Type-C Charging Cable *1

Packaging Box *1

Product Features Introduction

This product is a rechargeable multi-functional ultrasonic beauty device, featuring a combination of four functions: HIFU Ultrasonic, WARM Heating, EMS Microcurrent, and PULSE Pulse.

HIFU Ultrasonic Function: This function targets the SMAS layer of the human body with ultrasonic focus points, restructuring the collagen within the SMAS tissue without damaging normal tissue. It restores skin elasticity, smoothens fine lines, and enhances facial contours. (Indicator light is blue)

WARM Heating Function: The thermal energy emits micro-waves at a frequency of one million times per second, deeply warming the skin, promoting the breakdown of subcutaneous fat, and rapidly opening pores to facilitate skin cleansing. (Indicator light is red)

EMS Microcurrent Function: Various modes massage the facial skin, strengthening skin lift and firmness. (Indicator light is green)

PULSE Pulse Function: Mild heat opens the pores, enhancing voltage circulation deep within the skin and revitalizing collagen, thereby improving skin elasticity. (Indicator light is purple)

Operating Instructions

Step	Illustration	Operating Instructions
Power On		Long press the M button for three seconds to power on the device. If the display indicates a low battery, please fully charge the device before use.
Function Selection		Upon powering up, the display defaults to the HIFU function. Press the L key to cycle through HIFU, WARM, EMS, and PULSE functions. Confirm the desired function by pressing the M button to enter the working state. In this state, press the L key to adjust the energy level from 1 to 6. To switch functions during operation, press the M button to exit the current state, then use the L key to select a new function. Each session automatically begins a 10-minute countdown and exits the working state upon completion.
Care Areas		1. Whole Face Care: Divide the face into four sections, using the illustration provided in the care procedure guidelines. 2. Massage Head Usage: Keep the massage head moving slowly on the facial skin; avoid lingering in one spot for more than three seconds. 3. Gel Use: Each mode must be used with a gel. If the gel dries during use, reapply as necessary.
Shutting Down		1. Automatic Shutdown: The device will automatically turn off after 10 minutes of operation without any button interaction. 2. Manual Shutdown: While the device is on, long press the M button for three seconds; the display will turn off and the device will shut down immediately.
Cleaning and Maintenance		1. Cleaning: Wipe the device with a damp, wrung-out soft cloth. Avoid using alkaline solutions, alcohol, or volatile oils as these can damage the device's exterior and functionality. 2. Storage: After use, organize the device and its accessories back into the packaging box for next use. 3. Safe Placement: Store the device in a safe location and ensure no heavy objects are placed on it. 4. Battery Care: Do not use the HIFU and WARM functions when the battery is below 40%. Ensure the battery is fully charged before use. 5. Long-Term Storage: If the device is not used for an extended period, charge it before turning it on again.

Order of Care

Weekly Intensive Case

<HIFU> Mode: 10 minutes + < WARM > Mode: 10 minutes = Total Duration : 20minutes

①<HIFU> Mode: Apply the gel to the care area and activate <HIFU> mode.

② <WARM> Mode: Apply the gel solution to the care area and activate <WARM> mode.

※ Do not remain on the same spot for more than 3 seconds to avoid skin irritation.

※ Either mode can be used independently for a comprehensive 20-minute session.

	1 Glide from the chin along the jawline to the back of the ears, then proceed downwards from the neck to the shoulders.
	2 From the chin tip, glide along the facial contours to the temples, then from the nostrils across the cheeks to the temples. Start from the inner corners of the eyes, gliding along the eyes to the temples.
	3 Lift from the eyebrows to the hairline, moving in an upward direction.

★ Perform each of the three steps 2-3 times.

★ Maintain a slow motion and avoid stationary placement on any area.

★ Reapply gel or lotion if it dries during case.

Order of Care

Daily Care Routine

◆ <EMS> Mode: 10 minutes + <PULSE> Mode: 10 minutes

① Apply gel to the care area before starting.


Increase gel usage if skin feels dry during case. Insufficient gel may cause irritation and reduce case effectiveness.

※ Ensure the case head is in full contact with the skin.

※ Excessive case can stress the skin. Limit use to prevent multiple cases within a single day.

※ Clean any dirt off the case head with a soft cloth or tissue to maintain smooth operation.

※ Avoid using the device when feeling unwell.

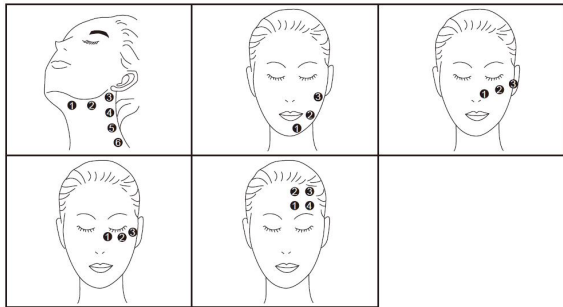★ Perform each step shown in Figure 1 2-3 times.

★ Place the device on each targeted skin area for about 3 seconds as illustrated.

※ Follow the illustrated case guidelines.

★ Apply more gel if it dries during the session.

Warnings

Do not use on mucous membranes, upper eyelids, eyeballs, areas around the ears, or the thyroid region.

Product Parameters

Product Name	HIFU Beauty Device	Model No.	HC096
Product Dimensions	180 x 60 x 54 mm	Package Size	225 x 108 x 90 mm
Product Net Weight	180 g	Packed Weight	390 g
Charging Voltage	DC 5V	Battery Voltage	7.4V
Battery Capacity	1500 mAh	Charging Current	500 mA
Power Rating	20W	Set Working Duration	10 minutes

Troubleshooting

Issue	Cause	Solution
Unable to Power On	Low battery	Please charge the device
Shuts Down During Use	Low battery	Please charge the device
Unresponsive During Use	Low battery	Please charge the device
Weak or No Current	Dry skin, poor conductivity	Apply a suitable amount of conductive beauty gel
Charging Failure	Incorrect charger voltage or current	Use the charger specified in the user manual
No Response When Charging	Faulty charger or charging cable	Replace the charger or charging cable

Warranty Card

Customer Name		Purchase Date	
Product Name		Customer Phone	
Customer Address			
Reason for Repair			
Maintenance Records			

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