
i3 Remeji

i3 Home-use Thermage

Instruction Manua

Before using the product ,please read this manual carefully and keep it properly.

This User Manual helps you understand the relevant functions and operation methods of "i3" more quickly.

Before using the product ,please read this manual carefully and keep it properly.

Please pay more attention to "Safe Use Instructions and Maintenance".



i3

All contents on this material have been carefully verified. In case of any printing errors or misunderstandings of the content, our company the right to modify the product design and specifications. Any changes to the product's appearance and color shall prevail over the actual product. No further notice will be given that time; Our company reserves the final interpretation of this "User Manual".

Note: This product is not a medical device category and does not have the function of directly used for medical beauty purposes

CONTENTS



01 Product Introductio

02 Product constitute

03 Product Accessories

04 Product show

05 product overview

06 Packing list

07 Process Intruction


08 Gear design

09 Instructions for safe use
and maintenance

10 Troubleshooting

11 After-Sale Warranty

PRODUCT INTRODUCTION

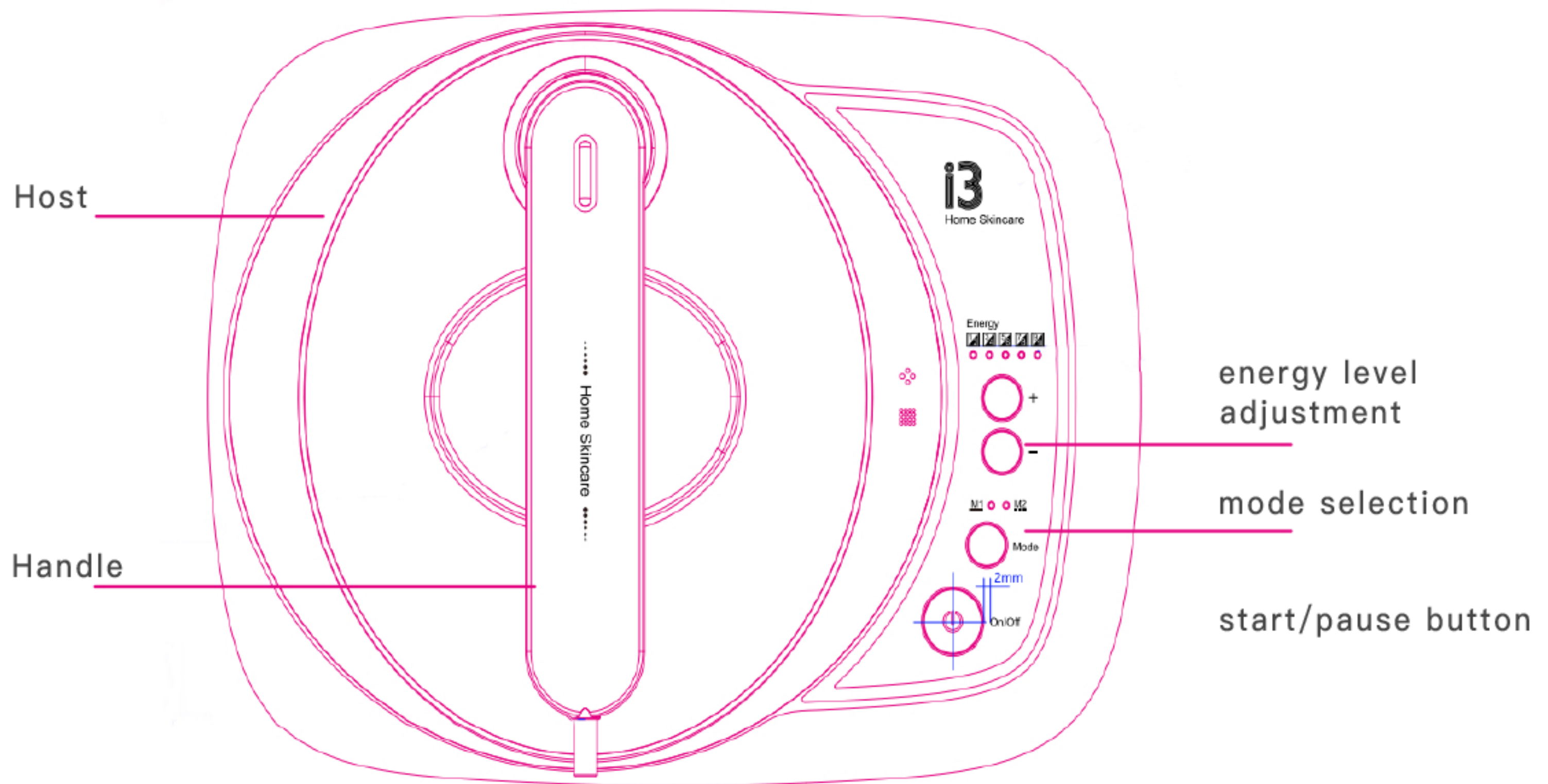

PRODUCT COMPOSITION

Main components	
Host	Power cord
Handle	4-pole array display
4-pole array	energy level adjustment
64-pole array	64-pole array display
mode selection	power cord interface
power switch	start/pause button

PRODUCT ACCESSORIES



Handle



64-pole array



4-pole array

PRODUCT DISPLAY



PRODUCT OVERVIEW

Basic parameter

Product model : i3

Product name : i3 Remeji

N. W. : 1.8 KG

Machine size : 32*24*7.8cm

Rated frequency : 50~60 HZ

Nominal voltage : 110V~220V available

Rated power : $\geq 80W$

Control output : The host button
and handle button
control adjustments

Language : English (languages can be customized)

Function : Skin Tightening and Lifting

PRODUCT LIST

Description	QUANTITY
Host	x1
Handle	x1
4-pole array	x1
64-pole array	x1
Instruction manual (including warranty card)	x1

PROCESS INSTRUCTION

- (1) Before use, thoroughly clean your face with a facial cleanser to remove surface dirt and pat dry.
- (2) Apply a thin layer of thermal insulation gel or essence cream on the area to be treated to reduce friction between the skin and the applicator head, enhancing comfort during use.
- (3) Select the 64-pin array head or the 4-pin array head based on the treatment area and treatment range. Align the array head the square and round interfaces of the handle and insert it. To change the array head, press the white button on the top of the handle to eject the array head. For delicate areas around the eyes, choose the 4-pin array head. For the face, neck, and other areas, it is recommended to choose the 64-pin array head.
- (4) Plug in the power cord, press the power switch, and the device will default to CM continuous energy output mode, with the energy set at level 1. It is recommended to use the default settings for the first operation of the i3, and then adjust the levels appropriately based on tolerance in later.

CM: Continuous power output mode (commonly used mode)

IM : Banded energy output mode (recommended for beginners or those unsure about their skin's heat tolerance)

注 : IM The band-like energy output mode is set with a human body induction system, and energy output only occurs after the array head contacts the skin following the of the handle

(5) Operating method

64 poles:

It is recommended to use a dot matrix stamping technique, with the probe lingering on areas that need care, such as nasolabial folds, jawline, and the "apple" area, etc., Each stamp should linger for 3–5 seconds before moving to areas for stamping care ; this process needs to be repeated. For areas that have been overheated, the stamping time can be appropriately reduced. Avoid staying the same position for too long to prevent overheating.

4 poles:

It is recommended to combine the dot matrix stamping technique with the sliding lift technique, such as stamping and dwelling on lower eyelid, then sliding and lifting from the lower eyelid to the temple. The eyes are more sensitive, the dwelling time is recommended to be 1–2 seconds, and the operation should avoid the upper eyelid and eyeball area.

(6) Operation time

The full face treatment time is recommended to be 10–15 minutes, and the eye area treatment time is recommended to be 5–8 minutes; the same area should not be treated more than three times a week.

- (7) After the procedure, the facial gel residue should be cleaned up in time to prevent pore blockage.
- (8) After the operation, turn off the power, unplug the power cord, clean the gel residue on the print head, disinfect it with an swab, and then dry it with a clean towel;; Wipe the handle and body with a clean soft cloth, and avoid using chemical cleaners or water for cleaning.

GEAR DESIGN

- (1) This product is designed with ten levels. The higher the level, the higher the temperature rise, and the greater the stimulation to the collagen fascia. Therefore, it is recommended for users: For daily care and skin without obvious problems, start more often from level 1, and adjust to a higher level according to the overall skin temperature during operation.
- (2) The heat should be set at a warm level, and you should not pursue higher settings for better results. For daily care, the setting is recommended be between 1–5. For those with obvious skin problems and experience using beauty devices, you can adjust the setting higher according to your own comfort level., the operation speed needs to be faster, and the stamping should be kept within 2 seconds to avoid excessive energy concentration and skin burns
- (3) Any heat setting recommendations should be chosen based on your skin's tolerance to prevent the risk of burns.。

Before using the instrument, please read the following contents carefully. If you do not strictly abide by the usage rules and precautions, the desired effect may not be achieved. If there is no progress in the place you want to change, or if you have any questions about this instrument, please consult your instructor as soon as possible. The results of use may be related to your age, living habits (exercise, diet, etc.) and attitude towards life.

- (1) Do not use this instrument on the part where jewelry such as jewelry or earrings are worn.
- (2) People with metal objects in their bodies (gold wire implants or steel nails, pacemakers, etc.) should not use this instrument.
- (3) People under the age of 18 are not recommended to use this instrument.
- (4) Do not use this instrument during female menstrual period.
- (5) Cardiac patients should not use physical therapy equipment.
- (6) If you have fever, infection, inflammation, skin diseases (eczema, herpes), warts, unscarred or newly scarred wounds, serious disorders or blood clotting problems, the use of this device is prohibited.
- (7) Do not use this device if you have phlebitis, just had venous surgery (or venous sclerosis) within two months, have severe Varicose veins or limb numbness. If in doubt, please consult your instructor first.
- (8) Pregnant women are prohibited from using this device.
- (9) Unless instructed by a physician, this instrument should not be used on the ear, inside the mouth, mucous membranes, or nipples.

- (10) Patients with epilepsy are forbidden to use this device. People with mental illness are prohibited from using this device.
- (11) If you take medicines such as anti-inflammatory, anticoagulant, antibiotic, antihistamine, etc. for a long time , or your skin will be bruised due to slight, Please consult your.physician before using this device.
- (12) After receiving any plastic surgery, you can use this instrument, but you must follow the doctor's instructions.
- (13) This instrument is only for professionals to operate.

warn

- 1 The operating head should always be kept clean.
- 2 Please do not use alcohol or acetone and other flammable substances when treating the skin before nursing.
If necessary, please use soapy water or medical distilled water.
- 3 If alcohol is used to clean and care for the bald head or to disinfect any parts of it, it must be soaked before working.

Electric shock protection: Class 1 \ Class BF

Anti-corrosion liquid: Do not use this machine in an environment where there is a mixture of flammable anesthetics and air, or nitric oxide.

:Operation mode :continuou

Precautions

- > Confirm that the operating voltage is 200~240V
- > Make sure the power cord is properly connected

- ① During the i3 Thermage treatment process, you may feel your skin heating up and turning red, and some people may experience pain. are all normal phenomena.
- ② This product must be used in conjunction with gel, cream, and other media. Using it with other essences, water, freezed powder, or using it while wearing makeup may cause the skin to produce a stress response.
- ③ Users should follow the recommended time and techniques in the manual during. They should not extend the usage time or the duration of the stamp stay to enhance the effect, as this may cause severe skin heating or even burns.

- ④ ant women, nursing mothers, those with skin wounds, inflammation, or allergies should not use this product.
- ⑤ During operation, no metal objects should be worn, those with metal implants in their bodies are prohibited from using it;
- ⑥ The energy level should be set at a warm and comfortable level, and should not be set high;
- ⑦ Absolutely avoid letting children touch the internal components of the machine through the small holes with metal objects to prevent the risk of electric shock

Custody:

- ⑧ This product is a home beauty device and cannot replace professional medical aesthetic treatments. If you have serious skin issues, please seek the help of a professional.
- ⑨ If you experience any discomfort during use, please stop immediately and consult a professional.
- ⑩ Do not place the product in a damp, high-temperature, direct sunlight environment, as it may affect the usage effect and safety.

FAULT HANDLING

1. The device fails to respond upon startup

- (1) Confirm whether the power plug is securely inserted into the socket.
- (2) Check if the socket is electrified, you can test the with other electrical appliances to see if it is functioning normally.
- (3) Inspect the power cord for any damage or breaks.
- (4) Try to restart device to see if it can return to normal display.
- (5) Gently shake the beauty instrument to check for any loose parts.

2. The handle has no heat output and no energy output.

- (1) Check if the connections between the handle, the connecting tube, and the main body are loose.
- (2) Set the energy the highest level and see if there is any warmth.
- (3) Check if the handle and the dot matrix head are fully inserted.
- (4) Increase the amount gel medium, set the energy level to the highest, and check if there is any warmth.
- (5) Check all the dot matrix heads and switch to different and levels to see if there is any warmth.

3. During operation, the machine suddenly stops.

(1) Check if the power cord is loose;

(2) Check if the probe is loose;

(3) Insufficient voltage.

If the listed failure phenomenon does not appear, please contact the after-sales service personnel to provide you with guidance.

(In order to allow you to use i3 Remeji more safely and efficiently ,Please accept the guidance and training of professionals first)

Routine maintenance:

- 1、 Instrument installation and use of the environment requirements clean, low dust, relative humidity $\leq 85\%$, no corrosive particles, conductive dust, etc, to avoid strong magnetic interference;;
- 2、 Choose a well-ventilated location, and do not place near heat sources such as radiators, thermal ventilation pipes or places where mechanical vibration will occur;
- 3、 Flammable substances are forbidden to be used for cleaning and disinfection of the operating head and operating parts during operation;
- 4、 Please keep all parts of the handle outlet clean to ensure the accuracy of the output energy;
- 5、 After use, the residual gel on the dot matrix head must be cleaned, disinfected with alcohol swabs, and then dried with a clean towel For the handle and body, please use a clean, soft cloth for wiping. Avoid using chemical cleaning agents or water for direct cleaning.
- 6、 Do not use gasoline, wine, diluent and other chemical agents treated with cloth to wipe the surface of the case;
- 7、 In order to ensure the normal use of the instrument, the purchase or replacement of the instrument accessories, please contact the manufacturer or agent directly.
- 8、 This machine is a precision instrument, strong vibration and impact should be avoided during movement and transportation, and inversion is strictly prohibited to prevent damage to the system and electronic components in the machine;

AFTER-SALE WARRANTY

The brand side provides long-term after-sales service for the instrument, life-long use guidance service, one-year warranty for the main unit, and three-month warranty for the handle:

- ① Host maintenance includes: host functional components, LCD screen failures caused by non-human factors;
- ② Handle maintenance includes: function operation handle, handle connection line failure caused by non-human factors.

The warranty period begins on the day the instrument is activated. After repair, if the same fault occurs again on the same part, it will be repaired free of charge within 30 days;The parts repaired and replaced during the warranty period belong to the company, and the parts repaired and replaced outside the warranty period belong to the user;The additional services and additional non-original configuration gifts promised by the agent shall be served by the original agent.

The following conditions do not meet the warranty conditions, and paid repairs will be provided:

- (1) Malfunctions caused by the user's negligence and improper operation;
- (2) Faults caused by natural disasters (fire, earthquake, flood, lightning strike, etc.);
- (3) Failures caused after the instrument warranty period expires;
- (4) Natural discoloration of relevant parts caused by heat generation in the characteristics of instrument components;
- (5) Free repair and replacement will not be provided for failures caused by disassembly and modification of the instrument by the user.

Warm reminder: In the event of a failure due to the above circumstances, the brand will still provide after-sales maintenance services, and reasonable maintenance will be charged for instrument or component failures Fees, cost will be charged if parts need to be replaced. The logistics costs incurred due to after-sales problems shall be borne by the user.

Warranty content

1 Parts of the main body of the instrument itself

2 function handles