

# Upgraded 4-in-1 ultrasonic knife

## Multifunction Beauty Instrument



## Device



## Product Functions

The upgraded 4-in-1 ultrasonic knife beauty instrument is composed of HIFU ultrasound /RF RF /EMS microelectricity /PULSE PULSE.

**HIFU ultrasonic function:** Ultrasonic focus on the SMAS layer of the human body, the collagen in SMAS tissue reorganization without damage to normal tissue, restore skin elasticity, fade fine lines, improve facial contour.

**RF RF function:** RF can send out 1 million microwaves per second, which can deeply warm the skin, promote the decomposition of subcutaneous fat, quickly open pores and help skin clean.

**EMS MICROelectric function:** MASSAGE facial skin in various modes, strengthen the skin lift and tighten.

**PULSE:** The micro-heat can open pores and promote the voltage circulation deep in the skin, make collagen more active and improve skin elasticity.

## Operation Method

- ① Insert the power adapter into the power socket, and the DC head into the instrument tail jack;
- ② Long press the power button for 3 seconds to enter the main interface.
- ③ Press L to select a mode. The HIFU/RF/EMS/PULSE can rotate downward to select a mode
- ④ Press the M key to enter the working state (20 minutes timing), when the mode is over or halfway can press the M key to exit the working state;
- ⑤ Press L key to select the energy size under working condition (1-5 cycle selection);
- ⑥ In the working state, press the switch button to close or open the head LED lights;
- ⑦ Long press the switch button for 3 seconds to shut down the machine.
- ⑧ After use, please pull out the power adapter, then pull out the power DC head from the instrument, and wipe the metal probe clean with a wet paper towel.




## Special treatment once a week

<HIFU> mode: 10 minutes +<RF> mode: 10 minute= Total 20 minutes

- ①Apply the gel on the part to be treated and perform a treatment in the <HIFU> mode
- ②Apply the beauty essence on the part to be treated and perform a treatment in the <RF> mode.

※Attention: It is not allow to stay on skin over 3s.

※Single Mode treament is allowed.

<div style="border: 1px solid black; padding: 5px; display: inline-block;">1</div> 	<p>Move along the face line from below he jaw to the rear side of the ear. Then , move from the parotid gland to the shoulder</p>
<div style="border: 1px solid black; padding: 5px; display: inline-block;">2</div> 	<p>Move along with the face line from the end of the jaw to the temple . Next , move along with the cheek from the nostril to the temple . At the end , move along below the eye from the inner corner of the eyes to the temple</p>
<div style="border: 1px solid black; padding: 5px; display: inline-block;">3</div> 	<p>Move from the bottom to the top while lifting the eyebrows upward to the hairline</p>

★Perform each procedure as steps 1 to 3 for 2-3 times.

★Move constantly and slowly and do not stay on the same part over 3s.

★Add more gel when it dries.

## Daily Treatment

◆ <EMS >Mode= 10 minutes

①Apply the gel on the part of treatment.

※During treatment, apply more gel if you feel dryness.


if using little gel, you may feel an irritation on the skin or the effect of the treatment may not be exhibited.

※Place all elements closely to the part of the skin to be treated.

Excessive treatment burdens the skin, Never use this product several times a day.

※During use, in case of that product's head becomes rough by the adhesion of dirt etc., wipe the surface of the head using the soft cloth.

※Do not use this product when you do not feel well.

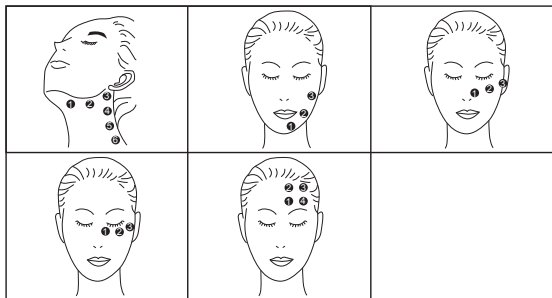

★Perform 2-3 times on each procedure.

★Place this product to the skin in order for approx. 3 seconds as shown in the illustration.

※The treatment procedure in the illustration is an indication.

★When gel becomes dry, add more.

### WARNING

Do not use one mucosa, the upper eyelids, the eyeballs, and around the ears and the Adam's apple

## FOR SAFE USE OF THIS PRODUCT

- Do not use the following person
  - People with heart disease, abnormal blood pressure, and malignant tumors
  - person who has a fever or an infectious or acute disease
  - People during dental treatment
  - Person who is suffering from a decline in immunity due to a deficiency
- disease such as AIDS or immunosuppressive drugs
- People with hemophilia, hemorrhagic disease, diabetes, facial neuralgia that have abnormal perception and pain
- Person with scar constitution, facial blackskin disease, photoallergy, pigmentation
- Chronic use of somatostatin drugs or telangiectasia caused by liver dysfunction
- People with allergic constitution, allergic dermatitis, skin diseases, sensitive muscles or fragile skin
- People who have a history of epilepsy or discomfort from flickering light
- People with head or eye sequelae caused by traffic accidents, etc
- Patients who have had eye surgery (consult your doctor before used laser in situ keratomileusis)
- Pregnant woman, physiology period, lactation period
- Young children, people unable to express themselves
- person who is drunk or hungover
- People who are taking sleep medications, extremely tired and unwell
- People who have had skin problems such as blotches and contact dermatitis caused by chemical fibers
- person who has or had dermatitis caused by cosmetics
- Recent chemical peels in skincare salons
- People who had skin problems after using beauty equipment
- Do not use in the following areas
  - Non-designated position
  - Mucous membranes surrounding the heart, eyeball, and upper eyelid, lips, and mouth
  - The site of a diseased or viral wart
  - The area of the skin that causes tingling and inflammation due to excessive sun exposure
  - There are wounds, eczema, dry spots, blisters, festering sores, and zits
  - recently been depilated or shaved

## PRECAUTIONS BEFORE USE

Prohibition	<ul style="list-style-type: none"> <li>● Do not use in conjunction with other cosmetic procedures.</li> <li>● Do not use for non-cosmetic purposes. Do not use in damaged or malfunctioning conditions</li> </ul> <p>This cosmetic instrument is a cosmetic nursing instrument and should not be used for therapeutic purposes or other purposes.</p> <ul style="list-style-type: none"> <li>● Do not peep or look directly at LED light</li> <li>● Do not use in the same place continuously. Otherwise it may damage the skin or cause skin problems such as infection, inflammation, bleeding, pigmentation, or worsen symptoms.</li> <li>● Do not use on pets.</li> </ul>
prohibition of use in a wet area	<ul style="list-style-type: none"> <li>● Do not use or keep in the bathroom or damp places</li> <li>● Do not place it where it is easy to touch water (washstand, etc.)</li> </ul>
prohibition of getting wet	<ul style="list-style-type: none"> <li>● Do not dip the device in water, or pour water over or sprinkle the device. (Otherwise it could lead to electric shock, fire or malfunction)</li> </ul>
prohibition of wet hands	<ul style="list-style-type: none"> <li>● Do not touch the device by wet hands. (Otherwise it will lead to electric shock and injury)</li> </ul>
prohibition	<ul style="list-style-type: none"> <li>● Do not use in the vicinity of diluent, volatile oil, deglazing water, spray and other combustible materials.</li> <li>● Do not use non-dedicated AC adapters</li> </ul> <p>The AC adaptor is exclusive to charge this product.</p> <ul style="list-style-type: none"> <li>● Don't connect the AC adaptor to other product, (Maybe it will cause fire, electric shock or failure).</li> </ul>
Unplug	<ul style="list-style-type: none"> <li>● Make sure to hold the end of the plug when plugging in or unplugging the adapter .Do not grab the cord to plug in or unplugging the adapter.</li> <li>● Do not connect the dedicated AC adapter to other devices</li> </ul>
prohibition	<ul style="list-style-type: none"> <li>● Do not damage the plug wire</li> <li>● Do not mutilate, bend, pull, twist, place heavy objects or clamp on it (Otherwise it may cause fire, electric shock or malfunction)</li> <li>● Please avoid overload of the socket (Otherwise it will cause heat or fire)</li> </ul>
Compulsion	<ul style="list-style-type: none"> <li>● Only applicable to household socket.</li> <li>● Please make sure the AC adapter and plug are tightly plugged.</li> </ul>



## CAUTIONS

Pull out the plug	<ul style="list-style-type: none"> <li>● Stop to use and turn off the power if the following condition occurs, remove the AC adapter from the socket</li> <li>● Earthquakes, thunder</li> <li>● This device or AC adapter, wiring smoke, abnormal heat, abnormal odor, deformation.</li> </ul>
Prohibition	<ul style="list-style-type: none"> <li>● This device may heat when using without any safety problems</li> <li>● Do not connect this device to metal or electrical products</li> <li>※ Mobile phones, televisions, etc.</li> </ul>
Prohibited remove	<ul style="list-style-type: none"> <li>● Do not disassemble, modify or repair the device by yourself and Please consult the After-service department when repairing.</li> </ul>
mandatory clause	<ul style="list-style-type: none"> <li>● Be sure to remove precious metals, glasses and contact lenses when using</li> <li>※ Precious metals: earrings, necklaces, watches, etc.</li> <li>● Be sure to start in 1 Level when you first use it or are not used to it yet</li> <li>● Do not press the head hard on the skin, do not stay in the same place for more than 3 seconds, otherwise it will cause skin problems.</li> <li>● When taking care of the skin around the eyes, please apply it to the eyes before use</li> <li>● If you feel abnormal skin or body pain while using it, please stop using it immediately and consult your doctor. Stop to use and consult a specialist</li> </ul>

※Please read the following items carefully before using

1	<p>Special care for each site = 1 time per week</p> <ul style="list-style-type: none"> <li>● Week 1-2 please nurse twice a week. The 3rd week after begins to inspect skin state and decide make care about every week 1.</li> <li>● Taking too much care is not good for your skin.</li> </ul> <p>Daily care should be performed at a frequency of 1 part per day or 2 days per day. Please use according to skin conditions.</p>
2	<ul style="list-style-type: none"> <li>● Please remove makeup and clean your face, remove makeup and sebum dirt, and apply special beauty lotion before facial care.</li> </ul>
3	<p>The state of the skin varies with the state of the body.</p> <ul style="list-style-type: none"> <li>● It depends on the condition of the skin. If the skin has rashes, redness, itching and other abnormal conditions, please stop using it and consult a specialist.</li> </ul>

## product specification

specification	Product Name	Four in one color screen ultrasonic knife
	Product Size	280*80*75MM
	Packing Size	
	Input Voltage	DC24V
	Working Current	600MA
	UF	2.4MHZ
	RF	3MHZ
	Electric Frequency	100HZ
	Impulse Frequency	200KHZ
	Timer	20minute
	Weight	480G

