



Safety Precautions

● Special attention

- Please assemble and use the product strictly according to the instruction manual
- Please read the instruction manual carefully and keep it properly for future reference.
- If you have any questions during use, please contact our customer service personnel.

 **Danger : personal injury due to incorrect operations**

 **Warning : danger due to incorrect operations**

 **Caution :personal safety or property loss due to incorrect operations**

DANGER:

- Under the normal working condition , do not insert any objects into the product , otherwise the product may harm the human body ;
- No data for use on children under 16 years old is currently available. Therefore, children under 16 years old are not allowed to touch or use it. Please keep away from children under 16 years old when using the product.
- The product should not be used by those (including children) with physical , sensory or mental impairments or lack of experience and knowledge of its use, unless supervised or directed by a person responsible for their safety.

WARNING

- * To avoid scald and failure when using , please regularly clean the dirt at the irradiation part ;
- * Please do not place the irradiation part near items that may easily cause demagnetization , otherwise it may cause failure and damage
- * Please do not use the product when the air outlet is blocked, otherwise it may cause smoke and fire ;
- * Please do not use the product on pets , otherwise it may injury ;
- * Please take a skin test before use , and confirm the skin condition after 2 hours . You can use the product if no abnormal conditions (redness or blisters) occurs in the irradiated area ;
- * After use , please wipe the hair and skin care products on the product to avoid failure and skin damage ;
- * Please set the irradiation intensity according to your own skin condition , and be sure to start with Soft Mode for the first time ;· It is normal that the skin is slightly hot and red after use ;
- * The product can not be left unattended after it is turned on , and please turn it off in time after use ;
- * If the operating environment temperature difference reaches more than 20 degrees (from low temperature to high temperature , or vice versa) , please turn it off and wait for about 2 hours before using it again ;
- * Please do not use the product near inflammable and explosive objects (including sprayers , gas stations , etc .) , which may cause fire , electric shock and so on ;
- * The product should be used at 0–45° C . Please confirm the room temperature before use ;

- * Please do not expose the product to sunlight or UV light for a long time , which will seriously damage the product ;
- * If the built-in fan does not work after the product is turned on (no humming sound) , please do not use it , otherwise it may cause smoke , fire , etc .
- * It is not recommended to use the product for a long time . please turn it off for 10 minutes after using it for 1 hour
- * If the casing , lens or power supply is damaged , please stop using the product immediately . Please do not repair the product without authorization . It must be repaired by professional from the manufacturer ' s maintenance department or similar department .

* If the operating environment temperature difference reaches more than 20 degrees (from low temperature to high temperature , or vice versa), please turn it off and wait for about 2 hours before using it again ;

* Please do not use the product near inflammable and explosive objects (including sprayers , gas stations , etc .), which may cause fire , electric shock and so on ;

* The product should be used at 0–45° C . Please confirm the room temperature before use ;

* Please do not expose the product to sunlight or UV light for a long time , which will seriously damage the product ;

* If the built-in fan does not work after the product is turned on (no humming sound), please do not use it , otherwise it may cause smoke , fire , etc .

* It is not recommended to use the product for a long time . please turn it off for 10 minutes after using it for 1 hour

* If the casing , lens or power supply is damaged , please stop using the product immediately . Please do not repair the product without authorization . It must be repaired by professional from the manufacturer ' s maintenance department or similar department .

● Unsuitable population and situation

* The product may cause discomfort or adverse effects on people with allergies , skin diseases , allergic dermatitis , sensitive skin , etc . If any adverse effects occur during use , please stop using the product immediately ;

* Please do not use the product on areas where the skin integrity is damaged (– e.g . , Burns , blisters , cuts , wounds , active skin diseases , or infections) ;

* If you are taking light – sensitive drugs , please seek medical advice before using the product , otherwise it may damage the skin ;

* The product is not suitable for pregnant women , lactating women and women in menstruation , with body scars , diabetes fragile skin and light sensitive skin ;

* People with severe heart disease , moderate to severe hypertension , malignant tumors , folliculitis , infectious skin disease , burned skin and anticoagulant disorders are prohibited from using the product ;

* Please do not use if you have abnormal collagen metabolism , including scars and poor wound healing ;

* Please do not use it if you have vascular diseases , such as varicose and hemorrhoids in the irradiation area ;

* Please do not use if your skin is sensitive to light and prone to rashes or allergies ;

* Please do not use it if your irradiated area is infected or you have herpes simplex , skin lesions and hematoma

* Please do not use it if you have immune response suppression diseases (including hiv infection , lupus erythematosus , porphyria , etc .) ;

* Patients with hemorrhage syndromes and those who are using isotretinoin , alpha hydroxyl acid (aha) , beta hydroxyl acid (bha) or azelaic acid are prohibited from using the product

* Please use this product carefully if you used chemical hair removal within 30 days or skin whitening products or cosmetics with vitamin a recently ;

* The product can ' t be used with essential oil products . Please ensure that the skin in the irradiation area is clean and dry without lotions , lubricants or other cosmetics on ;

* If the skin surface is not clean and dry , it may affect the user experience of the product and may even damage the product ;

* Those who are taking aspirin drugs , immunosuppressive drugs , isotretinoin and other drugs and recently took light sensitive drugs (such as tetracycline and sulfa) are prohibited from using the product ;

* Those injected with fillers and toxins recently are prohibited from using the product .

● Please do not use the product on any of the following body parts or any body parts with the following conditions

1/Oral cavity,genitals,mucosa,etc	5/Eyebrow,hair	9/Large vein
2/Throat knot	6/Ear	10/Tattoo,body painting
3/Eyeball,eyelid,around the eyes	7/Nipple,areola	
4/Head	8/Inside the navel	

(Note: This list is not exhaustive.If you are not sure whether you can use this product, please consult a doctor)

▲ NOTICE

1.The product is not intended for the treatment of diseases

2.Don' t expose yourself to the sun before using

3.Please use low level mode at your first use of the product and use the product on a small area on your arm to test your skin reaction

4.Make sure that the window is close to the skin before using the product

5.It is prohibited to use the product in dark light

6.To avoid failure , please do not use the product in direct sunlight

7.Please do not hit and shake the product violently

8.Do not unplug the product before turning it off , otherwise it may damage the product

Note : When you stop using the product , please place it with the air inlet upwards and the LCD screen downwards for better heat dissipation and cooling .

Product introduction



● About the product

*The product miniaturizes the medical beauty milk light technology. It heats the surface of the skin to proliferate and reorganize the collagen in the dermis, thereby achieving smooth and delicate skin. It can also stimulate the regeneration of subcutaneous collagen to achieve a firming and lifting effect.

*The product has an integrated lamp head that has a long life and does not need to be replaced, and is suitable for multiple body parts (such as face, torso, body, etc.), saving you expensive lamp head replacement costs, no need to purchase accessories frequently, and no hidden additional costs. In addition, it has a complete filtering system, which can effectively filter ineffective short wavelength and harmful light.

● PRODUCT STRUCTURE AND ACCESSORIES

PRODUCT COMPOSITION

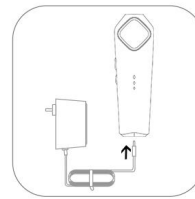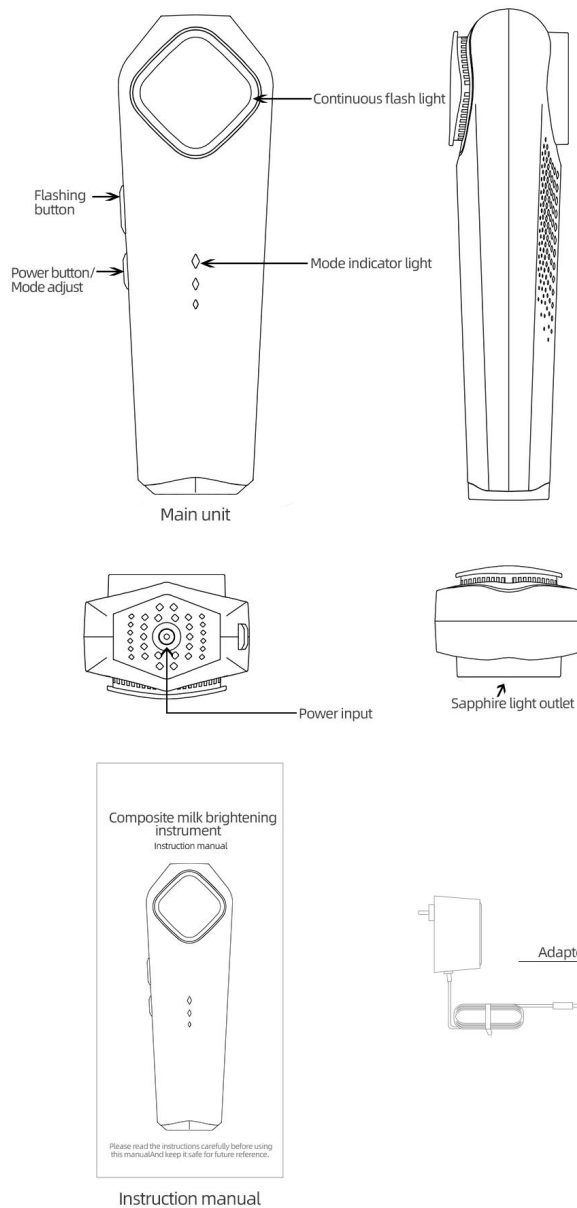

step1



Insert the plug into the power socket , after connecting the power cord to the product

*After finishing the operation , be sure to turn off the power before unplugging the power cord to prevent damage to the power supply

step2

Press and hold the power button on the side of the instrument to turn on the instrument

*After the instrument is turned on, the fan starts and the breathing light turns on

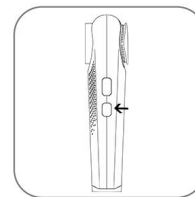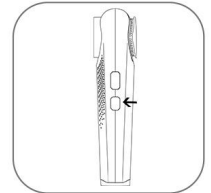

step3


Short press the switch key to select the appropriate gear

*There are three gears. It is recommended to use the first gear when using it for the first time

step4

Long press the flash button to switch between two modes: automatic mode and manual mode

*Automatic mode: After the instrument comes into contact with the skin, it will automatically perform care

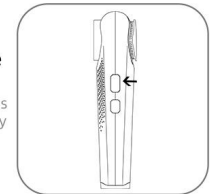

step5

Place the window of the instrument close to your skin and press the flash button to start the operation


*When the window is close to the skin, the indicator light will stay on

step6

After the care is completed, long press the power button to turn off the device

*After finishing the operation , be sure to turn off the power before unplugging the power cord to prevent damage to the power supply

*If it is turned on and not used, it will automatically shut down after 3 minutes.

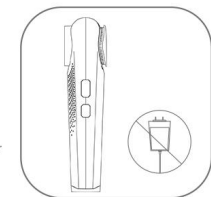

Preparation before use

● Before use



① Thoroughly clean the skin to ensure that it has no hair or oily substances (perfumes , lotions , deodorants , sunscreens , etc .) on and is completely dry .

② Make sure that there is no dirt , hair or other dirt at the irradiation port and that it is not damaged or has other abnormal conditions . If there is dirt , please clean it up to avoid failure , deterioration and skin problems .

③ Take a skin test before using the product and confirm the skin condition after 2 hours . You can use the product if no abnormal conditions (redness or blisters) occurs in the irradiated area ;

● Irradiation intensities

* The product has five modes with different energy intensities for different degrees of skin sensitivity

● About modes

* Smart continuous play mode

* Place the light outlet on the area to be treated and let it fully contact the skin . At this time , the sensor indicator light will light up , indicating that the lamp head fully contacts the skin . Long press the flashing button for 1.2 seconds and the continuous play indicator light will light up . Smart continuous play mode will engage . Please keep the light outlet contacting the skin . Flashes will be emitted at a certain interval . Please keep moving the product . Move it back and forth twice on the area to be treated to ensure that the area is covered evenly and avoid missing light .

* Tip : you can adjust your movement speed according to the flashing speed , so that the product can cover the area more evenly . The flashing frequency is higher in low level mode . Furthermore , the higher the level , the slower the light will be emitted .

● Single play mode

* Place the light outlet on the area to be treated and let it fully contact the skin .

The sensor indicator light will light up at this time , indicating that the lamp head fully contacts the skin . Short press the flashing button to start . Move the lamp head to the next area , wait for the white light of the sensor indicator light , and then press the flashing button again to repeat the above operation in order .

● Use cycle

Use 1-2 times a week

Care after use

● Precautions after use

- * It is recommended to apply a hydrating mask after facial treatment
- * It is recommended to apply body lotion after body care
- * If you feel discomfort , please stop using the product and Consult a doctor .

● Cleaning and maintenance

In order to ensure the best results and the longest service life of the product , please Clean the product before and after each use to avoid the long-term accumulation of hair and dust . If you do not clean it properly , the product will not realize its proper effects .

- * Before cleaning the product , you must turn it off , unplug it and let it cool down .
- * Use a soft , clean and dry cloth to wipe and clean the mirror of the light outlet .
- * Store the device in a dry and dust - free environment at -10-60° C .
- * Do not wash the product and accessories under a faucet or in a dishwasher .
- * Do not use steel balls , abrasive cleaners (toothpaste , etc .) or corrosive liquid (gasoline , motor oil , acetone , etc .) to clean the product .

Note : If the water can not clean the product well , you can use a few drops of high-concentration alcohol (lower than 96%) to clean the above parts .

● Environmental protection

When disposing of the product , please do not put it with general household waste . Please send it to a designated recycling center . This is good for the environment .

FAQ

This chapter outlines the problems you may encounter when using this product . If you cannot solve a problem with the following solutions , please contact the customer service center in your country / region .

- * When an abnormal situation occurs or a fault cannot be eliminated , please contact the after - sales customer service personnel in the corresponding language .
- * Please keep your purchase invoice and warranty card . You need to bring these certificates when the product fails .

Faults and solutions		
Problem	Possible reason	Solution
The product can't be turned on normally	The power supply is not connected	Check if the power supply is connected
The product can't emit light	The temperature of the main unit head and the body is too high The main unit head does not contact the skin closely	Wait for the device to cool down naturally before using it again; Make the main unit contact the skin closely until the light indicator light up
The working indicator light is not on	The product breaks down	Check if the power supply is working properly; Repair or return the product immediately
Odor appears when the product is flashing	The light filtration mirror is contaminated there is hair left on the skin surface	Wipe the light filtration mirror shave the hair and reduce light intensity
Unsatisfactory effect	You used a light intensity that is too low for you	Please choose a higher intensity next time
	You didn't press the flashing button when you moved to the next area	You should press the flashing button when you move to the next area and there should be some overlaps between different areas
	You did not set it after use Skin care	It is recommended to apply moisturizing mask to the face after use. Apply body lotion to body parts
	You used the product no frequently enough	the recommended irradiation regimen
	Your response to light exposure is slow	This product should be used continuously for at least 4 weeks
The product becomes hot and humid during use	Maybe the product has worked for a long time	Turn off the product and let it rest and fully cool down before using it again

Warranty regulations
<p>*The warranty period of this product may vary according to country in which it is sold , and we will try to extend it as for possible under the laws of that country .</p> <p>*The warranty period of this product , the power cord and accessories is 1 year from the date of purchase .</p> <p>*If the product has quality problems during the warranty period.We will repair or replace it free of charge . However , machine damage caused by improper use or improper operation is not covered by the warranty . We will replace the damaged part or the old device at the cost price . The old device will then be owned by us .</p> <p>*Please present your warranty card when you need the warranty service . The warranty date of the product is based on the date when you buy it . If there is no a date of purchase , it will be calculated by the date of production .</p> <p>In the following cases ,We may not provide the warranty service,refunding or changing it We is not responsible for reissuing the purchase invoice and warranty card for you . Please keep your purchase invoice and warranty card properly</p> <p>1.Product damage caused by unpredictable events such as fire , earthquake , flood , etc 2.Product failure and damage caused by your negligence , misuse or use under other abnormal conditions .</p> <p>3.Product damage caused by your unauthorized disassembly and improper storage (such as rodent damage or liquid penetration) .</p> <p>4.Product damage caused by the mandatory use of the product in an applicable environment , including but not limited to excessive workload .</p> <p>5.Product damage caused by other conditions that We cannot control .</p> <p>6.For damage associated with the function or use of this product (such as loss of your business benefits , business interruption , other losses caused by the use of the product , costs of replacing equipment , devices or services to use this product , your customer claims , etc .) , We assumes no responsibility .</p> <p>For the type and extent of accidents , We will bear corresponding responsibility according to relevant laws and regulations.</p>
<p>Quality assurance</p> <p>The company guarantees the product quality, reliable performance, safety and practical. In order to solve customers' worries,The company's implementation of quality responsibility, after-sales service, with this card to be handled.</p> <p>1. the product within one year after normal service, with this card in order to produce failure, free repair and replacement parts.</p> <p>2. The product warranty period is one year.</p> <p>Look out</p> <p>1. Abnormal use and artificial damage are not covered by warranty.</p> <p>2. If the customer disassembles the machine by himself, the warranty will be cancelled.</p> <p>Date of purchase: _____</p> <p>Where to buy: _____</p>