

Binding line



Home beauty device
instructions

Contents

1、 Product Overview and Parameters.....	01
2、 system interface.....	02
3、 Operation steps.....	03
4、 operating method.....	07
5、 Usage taboos.....	08
6、 precautions.....	09
7、 Instrument maintenance and upkeep.....	10
8、 warranty card.....	11

1、 Product Overview and Parameters

Welcome to use this home beauty device. Before use, please read this manual carefully and strictly follow the instructions in the manual to operate this beauty device.







working voltage : 28V	power: 28-56W Product	Product model: Z-2411
size: 21 * 4.2 * 4.6cm	Current: 1A	
Operation head: 4.5 heads per unit, 3.0 depth adjuster per unit, 2.0 depth adjuster per unit		
Applicable voltage of power supply: 220-240V		
Execution standards :		
GB 4706.1-2005 《Safety of household and similar electrical appliances - Part 1: General requirements》		
GB 4706.15-2008 《Safety of household and similar electrical appliances - Special requirements for skin and hair care appliances》		
GB/T 36419-2018 《Household and similar skin beauty devices》		

2、system interface




3、Operation steps

-  Firstly, clean the operating area thoroughly (if there is makeup, use makeup remover to remove it first, and then clean it with facial cleanser); 
-  Connect the power cord to the host. When the tpeyc interface of the power supply is plugged into the host, you will hear a clicking sound, indicating that the power socket has been connected in place. After the power plug is connected in place, the instrument will make a beeping sound of power on. Then press the circular button  with the power symbol in the middle of the instrument and hold it for 3 seconds to turn it on. Release the button immediately after turning on the instrument, and the display screen will light up.



Intensity is the energy intensity display, with instrument energy intensity ranging from 1 to 9 levels.
Shots used number display, (C0 unit is 10000. When the number of used posts at the bottom reaches 10000, the number will be displayed as C1, and 20000 posts will be displayed as C2)



Press the gear button  on the instrument to select the appropriate gear, and set the instrument gear to cycle setting, 1-9 cycles.

(For the first three uses of the instrument, it is recommended to perform testing operations from the lowest gear.) After selecting the gear, lightly press the button with the power symbol on the instrument (in the power on state, press the start button to release it without long pressing),



After the start button is pressed, the instrument is working and energy is being output. The instrument will make several beeping sounds,

After the sound ends, it represents that after one shot is played, the number of used shots displayed at the bottom of the screen will increase by +1;



4. Test energy intensity: First, pull at a constant speed from the lower right jaw position to the corner of the mouth. If there is no sensation, increase energy until there is a faint toothpick lightly touching the skin.

5. Sensitive individuals may experience slight heat, stabbing pain, or soreness during operation, which is normal. If the teeth are embedded with metal, they may feel a sudden increase in energy when operating the affected area. After passing the area, the energy returns to normal, which is also a normal phenomenon.

6. When operating near the nose, it is necessary to lower the energy level (adjust the energy level according to personal situation), because positions with more bones are more sensitive and can withstand lower energy (the energy level varies depending on the operating site, and needs to be adjusted according to one's own situation)

Normal skin: It is recommended to start using energy from a low level for the first time, and adjust the energy to the maximum personal tolerance based on the feeling during the operation; **Red bloodshot skin and sensitive skin**

with thin stratum corneum: Use low-end operation for the first time when using energy; Operate 2-3 times and then adjust the energy to the maximum range that the individual can tolerate based on the condition of their skin; Remember not to perform fixed point operations during operation. It is necessary to move slowly on the face, and operations are prohibited within the eye range and Adam's apple position.


7. Deep switching



1) There are two depth switches for the operation head, namely 2.0 and 3.0 switches. When the instrument does not use the switches, Depth of 4.5



2) Apply a layer of thick gel on the operating head, as shown in the figure, and then align the switch to be used with the socket on the operating head and directly cover it. After transparent part of the head film to bubbles or lack of gel, if there are no

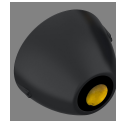

bubbles or lack of gel, You can start using it directly.

be careful:

When using the switcher, it is necessary to ensure that the gel of the round convex point inside the switcher is full. If the convex point covers the operating head, the gel between the operating head and the convex point is not enough and there is a gap, the energy transmission of the operating head will become weak or even no energy will be transmitted out. (It is recommended to watch the instructional video in detail before using the instrument)



8. During operation, the operating part needs to be coated with media (the media can be essence, stock solution, gel, aloe gel, lotion, etc.). During operation, the probe must be close to the skin, moving upward at a constant speed, and the media needs to be supplemented immediately after being absorbed or dried;

9. The normal operation time for the entire face is approximately 30-45 minutes; (Everyone's face size is different, so the operation duration is also different);

10. After using the instrument, clean the residual medium on the probe with a wet towel, and the gel in the switch can be directly cleaned with water.

11. After the operation, clean the face with clean water or warm water, immediately apply a water replenishing facial mask, and take it off after 10-15 minutes. If the water on the facial mask paper is absorbed, replenish the water in time; (It is not recommended to wash your face with hot water). After the operation in autumn and winter, your skin may be extremely dehydrated. Therefore, it is recommended to apply a pure moisturizing facial mask three days before the operation.

Attention: If you need to operate other parts, please consult a professional before proceeding.

4、operating method

neck



It is recommended to use the 3.0 and 2.0 operating heads for operation. When operating, it is recommended to start with the lowest energy. The probe should be lifted upwards along the neck from above the collarbone to the jawbone, then tightened towards the throat from the double chin position (not towards the Adam's apple), and then operated near the ear root lymph nodes from the double chin position; Remember to avoid the location of the Adam's apple: (During operation, the movement speed should not be too fast, and it can be operated back and forth 3-4 times, depending on the individual's relaxation situation)

face



It is recommended to use a mixture of 3.0 and 4.5 operating heads. When operating, it is recommended to start with the lowest energy. The probe should be pulled upwards from the jawbone to below the cheekbone, then pulled upwards from near the corner of the mouth towards the front of the ear. Finally, it should be pulled along the curved lower edge of the apple muscle from the nose to near the temple; (When operating, the movement speed should not be too fast, and can be operated back and forth 3-5 times, set according to personal relaxation)

eyes



It is recommended to use the 2.0 operating head for operation. When operating, it is recommended to start with the lowest energy. The probe should be pulled up from the frontal bone position to near the temple, and then pulled up from the back corner of the eye to near the hairline. **Remember not to operate into the upper eyelid, eye socket, and eyeball** (when operating, the movement speed should not be too fast, it can be operated back and forth 3-4 times, depending on personal situation)

forehead



It is recommended to use the 2.0 operating head for operation. When operating, it is recommended to start from the lowest point and pull the probe upwards from the center of the eyebrow to near the hairline, from the middle of the eyebrow to near the hairline, from the tail of the eyebrow to near the hairline, from bottom to top, and then from the forehead position to near the temple, and from the forehead to near the solar acupoint; It is recommended to avoid the eyebrow area (when operating, the movement speed should not be too fast, it can be operated back and forth 3-4 times, depending on personal situation)

• **Attention:** During operation, the probe should not be positioned in one location and repeated. The operation should be carried out at a constant speed. It is not recommended to adjust the energy to the highest level during initial operation. It is recommended to start from the lowest energy level. Adjust the energy level according to the condition of the operating part, and operate within the range of energy that the individual can accept; During operation, a slight heat sensation may be felt, and sensitive areas may experience a slight stinging sensation. After operation, there may be a feeling of soreness or numbness, which is normal;

- It is recommended to perform the first 5 operations on normal skin once every 5 days, and after 5 operations, it is recommended to use it once every 10 days, with 10 operations per cycle. It is recommended to perform the first 5 operations on sensitive skin once every 10 days, and after 5 operations, it is recommended to perform the operation once every 15 days. 10 operations are considered as one cycle, and the number of operation cycles should be set according to the degree of individual skin relaxation and sagging;
- After the operation, there may be varying degrees of detoxification metabolism on the face: acne, pimples, bumps, redness, itching, slight edema, soreness, etc. are all normal, and normal skin usually repairs automatically within 1-3 days. In some special cases, it may take a long time to repair. If the above situations occur, apply more water to the face. These are usually normal, mild, and temporary reactions.

5、Usage taboos

Individuals with the following symptoms should use this device with caution and seek the advice of a doctor or professional before use. The details are as follows:

1. Clean the operating area thoroughly before use to prevent dirt or debris. Before operation, remove any metal objects from the head to avoid unforeseen situations that may affect the effectiveness;
2. It is not recommended to use during pregnancy, postpartum period, or menstrual period;
3. It is not recommended to use the beauty device during the period of facial allergies. It is recommended to wait for a week after the allergy subsides before using it.
4. It is recommended to use the beauty device 20 days after facial injections such as hyaluronic acid or slimming injections;
5. It is recommended to avoid the area where prostheses and silicone are used for the operation site;
6. It is not recommended to use this area on the face where there are wounds, large areas of acne or pus;
7. Patients with serious disease, heart disease, diabetes, mental illness, etc. 8. Patients with serious skin disease, transfection disease, or hormone face should use caution;
9. It is not recommended to use it during pregnancy, but it can be used after postpartum confinement when the body has recovered. 10. It is not recommended to use it if the face has undergone gold wire implantation (metal wire cannot be absorbed);

6、precautions

- It is recommended to apply at least one water replenishing facial mask every day in the first three days after operation. If the conditions are sufficient, it is best to apply one water replenishing facial mask in the morning and evening of the first three days, apply at least 5-7 times a week of water replenishing facial mask, and drink no less than 2000cc (eight cups of water) every day;
- Eat less spicy food, smoke less, drink less, stay up late, eat more fruits and vegetables, and eat lighter food;
- It is not recommended to eat a large amount of soy sauce and other heavily pigmented foods within a week;
- It is not recommended to go to high-temperature saunas or other high-temperature environments within one week after operation, and do not expose to direct sunlight;
- Do not wash your face with hot water within three days (water that does not exceed body temperature is sufficient);
- No special care is required, normal hydration and sun protection are sufficient;
- During the operation, the operating area should be kept moist, and the probe should be kept tightly attached to the skin, otherwise energy will be concentrated on the surface of the skin;
- It is not recommended to operate on the eye socket, upper eyelid, temple or Adam's apple, nose bridge and other parts when using the face;
- When operating, it is recommended to gradually adjust the energy from small to large, and adjust it to the most suitable energy intensity for oneself. Do not directly turn on the maximum energy operation, as it may damage the skin (burn);
- The use of this product may cause discomfort or adverse effects on individuals with allergic constitution, skin diseases and allergic dermatitis, sensitive skin, etc. If adverse effects occur during use, the use of this device should be immediately stopped;
- Do not use this device in areas where skin integrity is compromised (such as burns, blisters, incisions, wounds, active skin diseases, and/or infections);
- Do not use this device on children.

Warning: Prolonged or repeated use of this device on a surface area may result in potential damage to the skin.

Warning: Do not bypass or cancel the inherent safety devices of the equipment, including cleaning and disinfection of the equipment or its components (if necessary)
The information.

7、Instrument maintenance and upkeep

1. After each use, please clean the instrument with clean water and store it properly.
2. Before using the instrument, it is necessary to use a plug with a grounding pin and ensure that the instrument power socket is truly grounded.
3. Ensure that the instrument voltage is compatible. If the local power supply voltage is unstable, we recommend that users add a power matching regulated power supply.
4. In order to ensure the effectiveness and normal service life of the instrument, please use the designated accessories provided or recommended by the original manufacturer uniformly.
5. The instrument should not be placed in damp places or near water sources, and should not be directly exposed to sunlight.
6. Do not place the instrument near strong heat sources, as this may affect the lifespan and normal use of the instrument.
7. Water, oil and other media must not penetrate into the interior of the instrument, and should not be hit, knocked or heavily dropped to avoid damaging the instrument.
8. Warranty period of the instrument manufacturer (within 15000 units). If the machine is damaged due to non-human factors during the warranty period, the manufacturer is responsible for free repair and replacement of parts (round-trip shipping costs need to be borne by the factory, and replacement and repair costs will be borne by the factory).
9. After the warranty period ends, the manufacturer provides lifelong maintenance for the instrument (replacement parts will only be charged at cost by the manufacturer, no maintenance fee will be charged).
10. It is prohibited to use items containing a large amount of alcohol for wiping and cleaning.

8、warranty card

Warranty Card	
Product name:	model:
Purchasing date:	Store name:
Address:	Tel:
Customer name:	wayyanty period:
Failure cause:	
Repair content:	
Warranty terms:	
1. If the product has quality problems caused by non-human damage, it will not be replaced unless there are no quality issues.	
2.If it is a software malfunction, the machine can be repaired within three days and will not be replaced.	
3. If the product has hardware malfunctions caused by non-human factors, our company will provide a one-year warranty from the date of purchase.	
Remark:	