



Bone
Hollywood International mode (Bone)



GOLDEN MUSHROOM HOT SKIN WAVE

GOLDEN MUSHROOM HOT SKIN WAVE MANUAL

KM-777 (Rose Gold)



Lily Singh



catalog

| | |
|---|---------|
| Product introduction | P.01-01 |
| Attention note | P.02-02 |
| Function introduction | P.02-02 |
| Machine introduction | P.03-03 |
| Handle detector introduction | P.04-04 |
| Control panel | P.05-05 |
| Installation operation | P.06-07 |
| Fault handing | P.08-08 |
| Machine maintain | P.09-09 |
| Machine parameter | P.10-10 |
| Parts list | P.11-11 |
| Progressive effect operation analysis | P.12-12 |
| Advantages Comparison of similar products | P.12-12 |
| Compare of before use and after use | P.12-12 |
| Machine certification | P.13-13 |
| Using safety and attention note | P.14-17 |

Preface

Golden mushroom hot skin wave is known as the most popular home beauty instrument now , it is suitable for girls aged 25 to 65.

In the early 1990s, American doctor Edward W.knowlton made great achievements on the study of RF acne removal and acne pit removal.Through research and found that the deep thermal effect of RF can kill more than 90% of Propionibacterium acnes in the dermis.By activating collagen and elastic fibers, the dermis tissue bitten by Propionibacterium acnes can be repaired, and the wrinkled and sagging skin can be restored to well state.

During the nursing experiment, the doctor's team found an unexpected receipt. After receiving the nursing care that changes the collagen at the part heated by energy, resulting in the contraction and tightening effect.

First. Product introduction


Introduction

Golden mushroom skin wave Lily Singh is a RF beauty tightening machine. When it applies energy to the skin tissue, the molecules in the tissue vibrate regularly, resulting in the generation of heat. This thermal effect will stimulate the renewal mechanism of collagen in the body, stimulate collagen tightening, promote the growth of new collagen, and achieve the effect of improving wrinkles, tightening and reshaping the contour.

Advantage

Safe noninvasive - Non invasive stimulation of deep collagen and elastic fiber proliferation
Deepest depth - Direct to deep fascia SMAS.
High accuracy - Focused dot matrix positioning, accurate operation without damaging surrounding tissues
One operation - Good nursing effect can be achieved without multiple operations
No recovery period - Not affect work and life after using
Lasting effect - Nursing effect lasts for a long time.



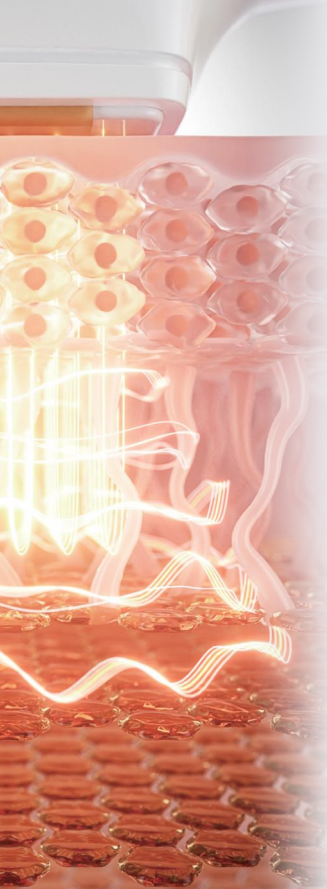

Second. Attention note

- 1.Keep skin and hair dean before operation. Do not make up before use. Make up traces need to be deaned.
- 2.It is necessary to lubricate the gel products before nursing, so as to ensure energy balance on both sides of the face.
- 3.After use, the skin will have temporary slight redness and swelling, which is normal.
- 4.Do not go to sauna, yoga and other high-temperature enviro-nment within one week after use. Exposure is not recommended.
- 5.Do not wash your face with hot water within three days (the w-ater temperature shall not exceed the body temperature); Pay at-tention to replenishment. It is recommended that the mask should be applied for at least three times a week.
- 6.Don't eat soy sauce and other foods with a lot of pigment within a week. Alcohol and spicy food are prohibited

Third.Function introduction

Golden mushroom skin wave Lily Singh is a skin care machine of high-safety and can't cause wounds.Its unique deep heating technology can stim-ulate the skin to regenerate new collagen, make your skin healthier and better, and improve the face contour. It is mainly used for overall skin tissue relaxation and sagging, excessive wrinkles, aging and rough skin, reshaping ace curve, and-d effectively improving the relaxation and sagging of tissues around eyes and neck.

- 1.Face: deep wrinkles,forehead lines, nasolabial groove lifting, lip wrinkles, eyebrow wrinkles, tightening neck skin, thin facial
- 2.Eyes: wrinkles around the eyes, crow's feet around the corners of the eyes, bags under the eyes, tightening the loose skin around the eyes



Wrinkles




relaxation



old age

Fourth. Machine introduction

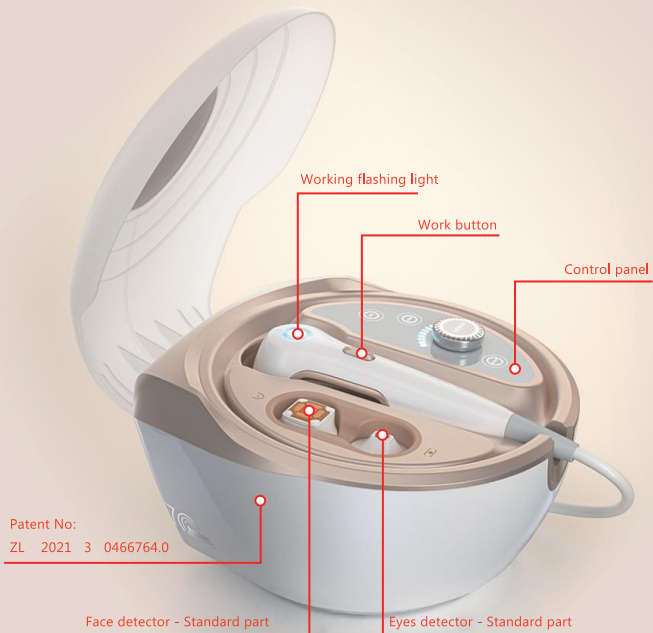


Fifth. Handle detector introduction



Handle Patent No: ZL 2021 3 0466777.8
Face detector Patent No: ZL 2021 3 0466776.3
Eyes detector Patent No: ZL 2021 3 0466743.9

Sixth. Control panel

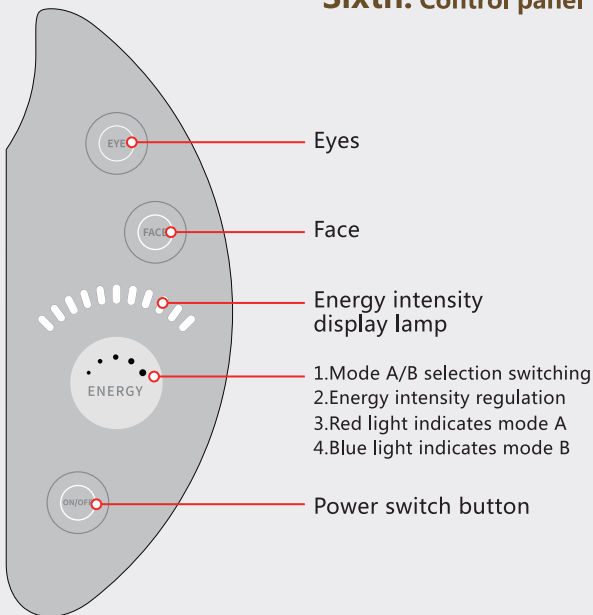

***Keying mode:** Press the button in the handle to start output,
release the button in the handle to stop output

Seventh. Installation operation

● Opening the box and checking

- Carefully take out the instrument and accessories from the packing box, and keep the packaging materials for transportation or storage in the future.
- Checking the accessories according to the accessories list and checking whether have any mechanical damage.

● Use steps



1 Cleaning

Before use, the caregiver is required to fully remove the makeup and ensure that the face and neck are clean.



2 Startup operation

Connect the power cord, gently press the ON/OFF button. On the host, Operation start, power off need unplug power supply.



3 Daub gel products to lubrication

Daub and make sure there is enough gel on the skin.



4 Adjust operating parameters

Start adjusting according to the energy adapted by the skin (it is recommended to adjust from a small energy).

Mode operation : After power on, rotate the ENERGY button to switch the A / B operation mode.
There are two operation modes;

(1) Mode A: The red light is on. Press the handle button to intermittently output energy.

(2) Mode B: The blue light is on, press the handle work button, and the energy output lasts for 1 second.
Press the button again, and the energy output stops.



5 Operation process

Light display above handle front end.



6 Cleaning skin

After use, clean the skin and use the replenishment mask for 10-15 minutes.

The operation detector is divided into: face operation detector and eye operation detector. Operate a single part (eye lines, legal lines) for 15 minutes, half of the face for 30 minutes and the whole face for 45 minutes. If the user feels that the heat is too large, reduce the energy. If there is slight thermal pain, it indicates that the pain is tolerable and the parts that need to care.

***Notes:** Since the operation mode is non-invasive, no recovery time is required, so daily activities can be resumed immediately

Seventh. Operating method

Operation production collocation: gel, Aloe Vera Gel



Neck

Stick to the skin from below the clavicle upward.
It is recommended to avoid the position of the Adam's apple.



Eyes

Use the fixed-point method to operate from the inner corner of the eye to the outer corner of the eye (it is not recommended to operate the part of the upper eyelid)



Forehead

From the middle of the eyebrows on both sides up to the hairline, and then from the position of the eyebrow peak, eyebrow ridge and the eyebrow tail to the hairline.



Cheek

Operate on the face according to the "shaped eyebrow ridge and the eyebrow tail to the hairline."

1) From the cheek, the operation of fixed point is carried out in three lines from the inside to the outside.

A Operate from the mandible to the earlobe **B** Operate from the corner of the mouth to the middle of the ear

C Operate from the nose to the temple

2) Then from the cheek, the fixed-point operation is carried out in three lines from bottom to top

A Operate from the mandible to the lower part of the frontal bone **B** Operate from the middle of the mandible to the lower part of the frontal bone

C Operate from the masseter muscle to the direction below the temple

3) Finally from the cheek, it is divided into three lines from inside to outside and operated by lifting and pulling

A Operate from the mandible to the earlobe **B** Operate from the corner of the mouth to the middle of the ear

C Operate from the nose to the temple

*Attention:

- 1, The operating head shall be close to the skin during operation.
- 2, The position with less fat such as cheekbones, frontal bones and eyes has obvious heat feeling. It is recommended to reduce energy.
- 3, When replacing the operating detector, push the pushing and twisting part, and the detector will pop up automatically. Do not pull it hard.
- 4, After each operation, the operating detector needs to be cleaned with warm water or normal saline.

Eighth. Fault handing

| Fault type | Fault cause reason judgment | Troubleshooting methods |
|--|--|--|
| After the power is turned on, the host has no indication and the operation button has no response. | <ol style="list-style-type: none"> 1.The power supply is not energized. 2.Damaged wires or poor connections. 3.Machine internal circuit failure. | Check the power supply to confirm power on or try to replace the power cord. |
| The working time indicator flashes (once every 1s). | <p>Flashing indication when the host cannot detect the operation head</p> <ol style="list-style-type: none"> 1.The operating head has been removed. 2.The operating head has poor contact. 3.The operating head is damaged. | Remove and insert the operating head repeatedly or replace the operating head. |
| In case of other faults, please contact your local dealer for repair or replacement | | |

***Attention :** The instrument is a precision electronic product. It is forbidden to use the instrument in a humid environment, such as bathhouse, bathroom and other humid places. If the instrument enters water, please stop using it immediately and return it to the factory for maintenance.



Ninth. Machine maintain

- 1.The instrument must use the original plug and ensure hat the instrument power socket is in good condition.
- 2.The instrument shall not be placed in a wet place or near the water source, nor shall it be directly exposed to the sun.
- 3.The instrument shall not be close to strong heat objects or strong radiation objects (disinfection cabinet, induction cooker, microwave oven, etc.), which will affect the normal use of the instrument.
- 4.The instrument must use the original qualified accessories provided by the manufacturer.
- 5.After operation, the operating head shall be clean, disinfect and properly keep.
- 6.When the instrument is not in use, please turn off the power supply, unplug the instrument and place various accessories of the instrument.

Tenth. Machine parameter


| Host parameters | | | |
|-----------------|-------------------------------|------------------|-------------------|
| Product | Golden mushroom hot skin wave | Frequency | 700-1060hpaeffe |
| Model | KM-777 | Working humidity | 10-90% |
| Electric source | AC 90V~240V 50HZ/60hzpov | Machine weight | 2.8KG |
| Power | 150W | Packing size | 367 x 323 x 235mm |

| Face detector parameters | | | |
|--------------------------|--------|--------------------|------------------|
| RF | | LED | |
| Power | 4J/cm2 | Wavelength | 650nm |
| active area | 3cm2 | Strength | 5000mcd |
| Frequency | 4MHZ | Energy output mode | Keyboard control |
| Active depth | 3mm | | |

| Eyes detector parameters | | | |
|--------------------------|--------|--------------------|------------------|
| RF | | LED | |
| Power | 4J/cm2 | Wavelength | 650nm |
| active area | 1.5cm2 | Strength | 5000mcd |
| Frequency | 4MHZ | Energy output mode | Keyboard control |
| Active depth | 3mm | | |

Eleventh. Parts list

| | Product | Quantity |
|---|-------------|----------|
|  | Machine | 1 |
|  | Handle | 1 |
|  | detector | 2 |
|  | Hangbag | 1 |
|  | Certificate | 1 |
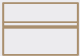
| | Product | Quantity |
|---|---------------|----------|
|  | Manual | 1 |
|  | Power cord | 1 |
|  | Packing box | 1 |
|  | Warranty card | 1 |

Twelfth. Progressive effect operation analysis

After use: Slightly lift + improve skin texture (obvious firmness)

First month: Slightly lift + improve skin texture (obvious firmness + Smooth skin)

Second month: Significant lifting + skin improvement (firm face + improved skin around eyes + obvious contour)

Thirteenth. Advantages Comparison of similar products

1. The golden mushroom hot skin wave can act at a deeper level and more accurately: RF has strong penetration and can accurately act on the dermis and 4.5mm SMAS layer. In terms of service temperature, RF can accurately focus on the dermis and SMAS layer, and the effect is more obvious.

2. The golden mushroom hot skin wave is safer to use.

3. Golden mushroom hot skin wave is non-invasive, mild and more comfortable, no irritation to the skin and slight pain.

Fourteenth. Compare of before use and after use



Sixteenth. Using safety and attention note

Professional tips:

- Before use, please read the manual carefully in order to use the instrument correctly.
- The precautions are for the convenience of the operator to use the instrument safely and correctly.
- In order to indicate the degree of harm or damage, the misuse in the precautions is divided into warning and attention.

In order that you can operate this instrument normally :

This manual is marked by diagrams and icons, and their representative meanings are as follows :

- After understanding the meaning of the logo, please read this manual
- After reading, please keep this manual in a place you can find at any time
- As the following contents are related to safety issues, please strictly abide by them



Warning

Possible injury or damage to the instrument

Meaning of icon symbol



Attention

There will be danger, warning and attention in the prompt, and the specific prompt contents are recorded in the icon. (the figure on the left is applicable to no specific and general hazards, warnings and notices)



Warning

There are prohibited behaviors in the prompt, and the specific prompt contents are recorded in the icon and below the icon. (the figure on the left is applicable to no specific and general prohibition notice)



Force + Tip

Prompt the behavior that needs to be enforced or instructed, and record the content of specific coercion or instruction in the icon (the left figure is applicable to no specific general coercion or instruction)

Taboo matters

Warning



This instrument cannot be used together with the following medical electronic instruments

- Pacemaker and other implantable medical electronic instruments
- ECG scanner and other installed medical electronic instruments
- Artificial heart lung life sustaining medical electronic instruments

Warning

This product is only suitable for beauty

Use for other purposes may cause accidents or obstacles

Don't use the product about the following people :

- Patients with heart disease and lupus erythematosus ● Light allergy, light syncope ● Patients who are taking antibiotics for acne (such as tetracycline and isotretinoin)
- Nervous and weak visceral function ● Menstrual period ● Suffering from acute disease ● Patients with atopic skin inflammation, high fever and fever, who have just received injection therapy
- Suffering from severe skin disease ● Diabetic ● Allergic constitution ● Fracture or metal joint ● Being ill or hospitalized ● Abnormal blood pressure ● Suffering from infectious diseases
- Hemophilia ● Suffering from tuberculosis ● Suffering from malignant tumor ● Suffering from facial neurosis ● During drinking or taking medicine ● Sharp fatigue ● Pregnant
- Telangiectasia caused by long-term use of steroids and liver dysfunction ● Suffering from legal infectious diseases ● Suffering from facial leukoplakia ● Children and adolescents

Be careful to cause accidents or adverse reactions

Please do not use this product in the following parts:

- Patients with skin diseases ● Pain and perceptual disorders ● Site with acute infectious disease ● Eyes and skin around eyes
 - Where metal, plastic and silica gel are implanted into the face ● Oral mucosa ● Acne, abscess and other inflamed parts
 - Parts with facial black spot ● Parts of lesions caused by itching, burning, physical stimulation, etc ● Wound
- Parts of skin allergy and inflammation caused by cosmetics ● Extremely dry skin ● Parts of skin irritation and inflammation caused by excessive sun exposure
 - The place where metal, plastic and silica gel are implanted in the body ● Anus and pudenda ● Head and hair

Be careful to cause skin problems

❶ Please do not carry out overload beauty care.

Be careful to cause accidents or other problems

❷ Please unplug the power plug immediately when the following abnormalities occur:

- When metal, water or foreign matters enter the instrument
 - The outer package is damaged due to falling, bumping and other reasons
 - When the product smokes, produces peculiar smell or noise.
- Be careful to cause fire or other problems

❸ Keep out of reach of children.

Be careful to cause accidents or other problems

❹ Plug in the power supply first, and then turn on the main switch of the instrument.

Avoid current backflow

❺ Do not use damaged plugs and loose sockets.

Be careful to cause fire, short circuit and electric shock

❻ Make sure the accessories are connected to the host.

Avoid socket damage caused by improper connection of accessories

❼ After use, be sure to pull the power plug out of the socket.

Be careful to cause fire or electric shock

Warning I



- Do not use or store this instrument in bathrooms and places with large moisture and humidity. Beware of internal water ingress, failure, fire or electric shock
- If there is water inside, please stop using it immediately. Be careful of electric shock
- Please do not plug and unplug the power plug with wet hands. Be careful of electric shock
- Please do not use wet hands to operate the instrument or load and unload the socket. Be careful of electric shock
- Please do not repair, disassemble or transform the instrument host and accessories without authorization. Beware of fire or scald caused by improper operation
- Please do not use this instrument in places with a lot of dust. Be careful to cause fire and instrument failure
- Please do not knock or puncture the operation panel of the host with sharp objects. Be careful to cause failure
- Please do not use or store this instrument in places with direct sunlight and near heating equipment such as stoves. Be careful to cause failure
- Please do not impact or drop the instrument violently or place heavy objects on the product. Be careful to cause failure
- Please do not use this instrument with other electrical products. Beware of misoperation
- Please do not use this instrument near computers and other precision instruments, televisions, mobile phones, telephones, etc. Beware of misoperation
- Please avoid storing the instrument in an environment that may produce strong electric waves and strong magnetic forces. Be careful to cause failure
- Please do not place containers containing water and other liquids on the instrument host. Be careful to cause fire, electric shock and failure
- Please do not use the spare parts of this instrument for other instruments. Be careful of scalding and accidents
- Please do not let children operate this instrument. Be careful to cause accidents or other problems
- Please do not use special power cord or other unexpected wires. Be careful to cause fire or failure
- Please do not come into direct contact with air containing salt and sulfur. Be careful to cause failure

Warning II



- Please do not use this instrument for other purposes and purposes. Be careful to cause accidents or other problems
- For the maintenance of the main machine and accessories of the instrument, please refer to the operation manual. Be careful to cause accidents or other problems
- Please stop using this instrument immediately when snowing. Be careful to cause abnormal output power
- Do not leave the customer during project care. Be careful to cause accidents
- Be sure to conduct initial inspection when purchasing. In case of abnormal conditions, check and repair immediately. Be careful to cause accidents or other problems
- Please check the power plug and socket regularly. Beware of fire or fire events
- Please regularly remove the dust covered on the power plug and socket to avoid fire or fire

Attention

Precautions for use

- The instrument can only be operated after knowing and mastering the correct use method of the instrument
- Cooperate with the company's designated skin care products for instrument care
- Please take care step by step according to the conclusion of the customer's skin test and the degree of skin needs. When the customer's skin and physique reach the tolerance limit, please adjust the mode in time
- Please confirm that the accessories and electrical specifications of each part of the instrument are in good condition before use
- If there is a sudden power failure during operation, the power supply must be cut off in time to prevent the instrument from being burned out by a sudden power call
- After nursing, make sure to cut off all power supplies related to the instrument

Attention

- When installing the instrument, please be careful not to drop the instrument during loading and unloading



Lily Singh

GOLDEN MUSHROOM
HOT SKIN WAVE

KM-777