


# Handheld beauty device instruction manual

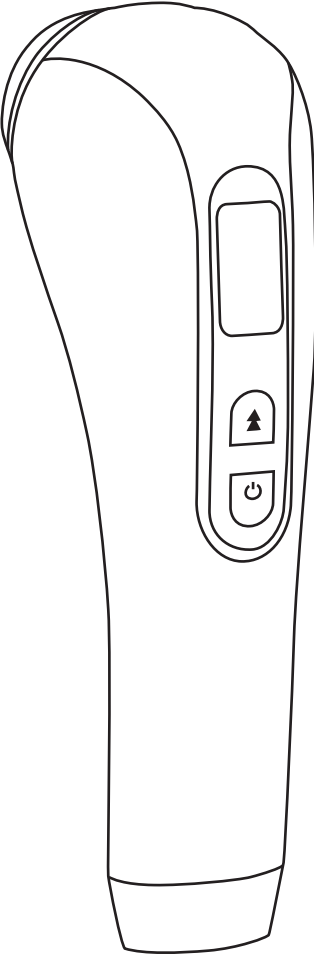

# Table of contents

---

01.Product Overview.....	02
02.Instrument parameters.....	02
03.Product diagram .....	02
04.Operation steps .....	03
05.Contraindications for use .....	05
06.Things to note .....	05
07.Instrument maintenance .....	05
08.Packing list.....	06
09.Warranty card .....	07

## 01.Product Overview

---

Welcome to use this beauty instrument. Before using this instrument, please read this manual in detail and strictly follow the instructions in the manual to operate this instrument.


## 02.Instrument parameters

---

Product Name	Handheld water drop beauty device		
Product model	H-0251	Power	96W
Power supply input voltage	100-240V	Frequency	3/10MHZ
Instrument working voltage	24V	Operation head	1
Current	4A	Instrument size	25*25*10cm




## 03.Product diagram

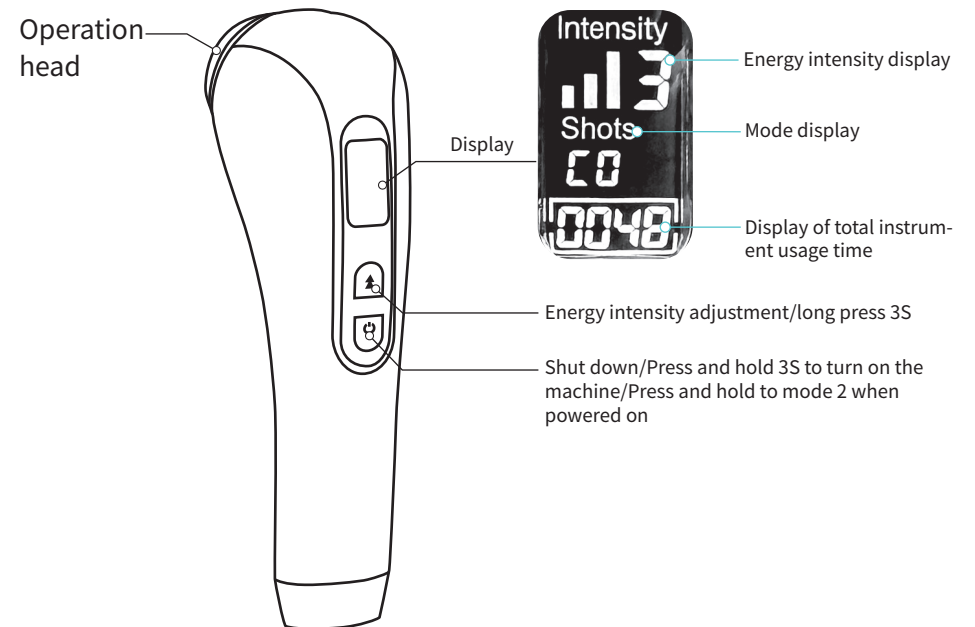





---



## 04.Operation steps





**Step 1:** First clean the operating area.

**Step 2:** Connect the power cord , note that the plug should be inserted into place .after the power supply is connected, insert the power supply tail into the power outlet at the tail of the instrument .after connecting to the power supply:

- **1. Power on:** Press and hold the power symbol button on the instrument  after 3 seconds, the instrument is turned on and the screen light is on.
- **2. Power on/pause:**Power symbol button  it is the start and pause button. When the power is on, press the button and the instrument will start. When the instrument is working, the instrument will pause.
- **3. Mode selection:** When power is turned on, press and hold the power symbol button  3 seconds, Shots on the display flashes, the current mode of the instrument is Mode 2 (Mode 2 is 3M sound wave output for 6 seconds, after 1 second pause, 10M sound wave output for 6 seconds, after 10M sound wave output, 1 second pause, and then 3M sound wave starts to output the cycle working mode).keep pressing the power symbol button in mode 2  for 3 seconds, the Shots on the display stop flashing and the instrument enters mode-operation mode.
- **4. Adjust gear:**Arrow button on the instrument , it is the gear adjustment button. When the power is on, press the button and the gear will change from 1-9 cyclically. Just select the gear you choose according to your needs.
- **5. Shut down:** If you need to shut down, press the button of the arrow , press and hold for five seconds, the instrument will turn off and the screen light will turn off.

### Step 3:Detailed operation explanation

● Apply special media on the operating area (the effect of removing acne is used with acne removal medium, and if it is used with acne pit removal or anti-allergicity, it is used with the corresponding media. The media can be gels, essences, raw liquids, creams and other skin care products). Apply the media evenly to the operating area. It is the limit that the white media is not visible. The media should be replenished immediately after the media is absorbed during the operation. During the operation, you need to slowly adjust the energy intensity according to your own feelings, gradually increase from low to high, and adjust the most comfortable energy intensity within the range you can accept to operate.

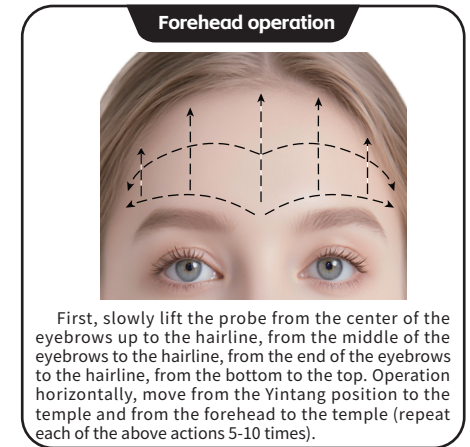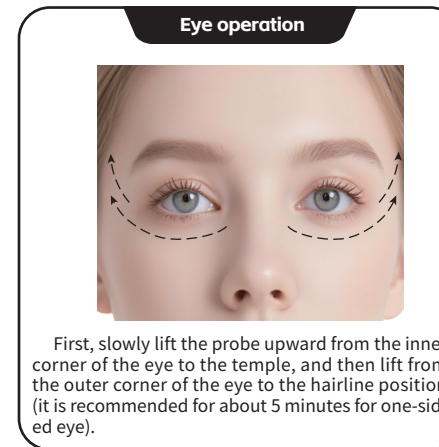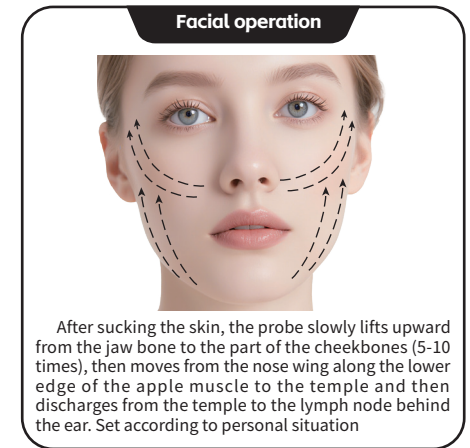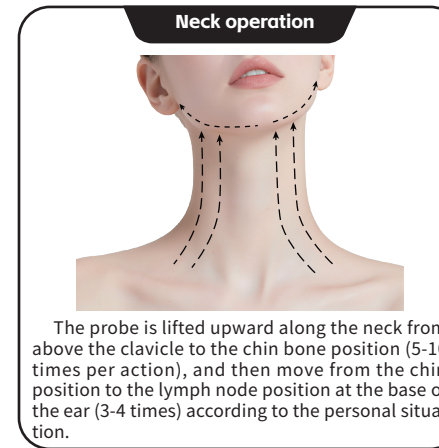

**Step 4:** After the operation is completed, clean the media in the operating area and apply a pure moisturizing mask or eye mask for 15 minutes.

**Step 5:**Use an alcohol cotton swab (without dripping as the standard) to wipe and disinfect the probe to prevent media from remaining on the probe and extend the service life of the machine.

### ● Operation time period:

It can be operated once a day, and the whole face can be operated for 15-20 minutes each time. For those with allergic skin, it is recommended to use it with your existing skin care products. When using allergic skin, you must operate at low energy for the first five times (not recommended for those who are allergic to metals).

**⚠ Notice:**

Everyone's skin has different skin abilities, and the energy intensity used during operation is also different. During the operation of this instrument, you must flexibly adjust the energy intensity of the instrument according to the actual feeling of your skin.

## 05. Use taboos

**Those with the following symptoms should be careful when using this instrument. Please consult the doctor or professional before use. The details are as follows:**

- 1 If the operating part is prosthetic, avoid this area when operating the instrument;
- 2 There have been injections of hyaluronic acid and collagen in the operating area within the past six months;
- 3 It is not recommended to operate if there are cardiac diseases and the device has a rhythm adjuster;
- 4 Do not put on makeup on the day of operation, or clean your face thoroughly before operation;
- 5 Avoid operations during women during menstruation;
- 6 It is not recommended to use this instrument for people allergic to metals.
- 7 Although there is no evidence to prove that this instrument will have an impact on patients with skin with tumors and active skin diseases, it is still not recommended for this population to use this instrument.
- 8 Although this instrument has not been directly proved that it is harmful to pregnant women, it is still not recommended for pregnant women to use it.

When using this beauty instrument, remember that energy must be used from the minimum and slowly adjust the energy up according to your personal feelings until your skin is most comfortable to operate! Do not turn on the energy directly to the maximum operation. Everyone's skin bears different energy, and the energy intensity selected during operation will also be different. When you operate, the instrument will directly turn on the maximum energy operation, which may cause damage to the skin!

## 06. Things to note after operation of the beauty instrument

- 1 Strengthen the moisturizing and sun protection of the operating area;
- 2 Please do not put on heavy makeup immediately after operation. It is

best to have more than 3 hours apart from the operation and makeup;  
3 It is recommended to apply a facial mask after each operation and drink no less than 2,000cc of water (eight cups of water) per day.

## 07. Instrument maintenance

- 1 Please clean the instrument with clean water and store it properly after using it.
- 2 Before using the instrument, you must use a pinned plug and ensure that the instrument power socket is truly grounded.
- 3 Ensure that the instrument voltage is adaptable. If the local power supply voltage is unstable, we recommend that users increase the power matching stable power supply.
- 4 In order to ensure the effect of the instrument and ensure normal service life, please use the designated accessories provided or recommended by the original manufacturer.
- 5 The instrument should not be placed in a humid place or near a water source, nor should the instrument be exposed directly to the sun.
- 6 Do not move the instrument close to a strong heat source, because this may affect the life and normal use of the instrument.
- 7 Water, oil and other media must not penetrate into the instrument, and must not be impacted, knocked or severely dropped to avoid damage to the instrument.
- 8 Please clean the probe in time and disconnect the power supply after use.
- 9 The warranty period is one year. During the warranty period, the machine is responsible for free repair and replacement of parts under normal use.
- 10 Provide lifelong maintenance and technical support.
- 11 It is prohibited to use alcohol-containing things for cleaning and wiping.

## 08. Packing list

Name	Host	Power cable	Power supply	Certificate of Conformity	Instruction	Warranty card	Test report
Image							
Quantity	1 unit	1 item	1	1 serving	1 serving	1 photo	1 photo

## 09.Warranty card

---

Maintenance record table		
Product name: Handheld water drop beauty device	Purchase Date	Warranty period of one year (from the date of purchase)
Dealership name:	Address:	
Customer name:	Telephone:	
Repair date	Cause of failure	Repair content