


3 in 1 Home Use Beauty Device User's Manual

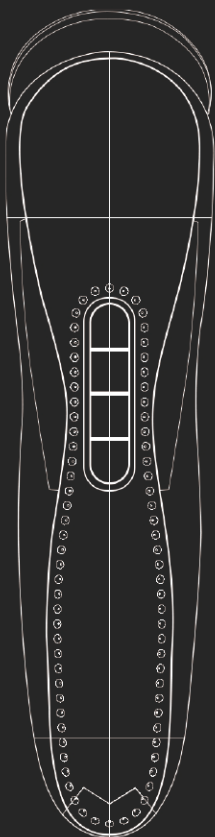

Model: TRI-ELF

*Please read this User's Manual carefully for safe and proper use.

TABLE OF CONTENTS

CONTRAINDICATION & SAFETY INFORMATION	2-7
Instructions	2
Safety precautions	2-5
Contraindications	6-7
WHAT IS THE TRI-ELF?	7-8
Definition	7
Technology	7
Story	8
WHAT'S INSIDE THE TRI-ELF BOX ?	8
GETTING TO KNOW YOUR TRI-ELF DEVICE	8-10
General guidance before use	8-9
Product structure	9-10
TREATMENT PROCEDURE	10-17
First time use	10
Basic notice	11
Mode & level adjustment	11
Treatment modes recommendations	12
Use & operation method	13-16
LED	16-17
CLEANING AND STORAGE	17
TECHNICAL SPECIFICATIONS	17
WARRANTY TERMS	18
Warranty terms	18
Warranty card	18




1.1 Instructions:

- ◆ Thank you for purchasing BCD product.
- ◆ Make sure to include this User's Manual if you are transferring ownership of the device to others for safety purpose.
- ◆ Make sure to store this User's Manual in an accessible location and fill in the information of purchasing date and place in the warranty card enclosed.




1.2 Safety precautions:

- ◆ Read all instructions and cautionary markings in the user's manual and on the device before use and always follow.
- ◆ Immediately discontinue use and arrange to have the device repaired if any abnormality occurs.
- ◆ Do not use the device in a defective state, as doing so could cause a physical harm to the user.
- ◆ The precautions listed in the User's Manual are very important. For your safety, please follow them carefully.
- ◆ Precautions are classified as DANGER, WARNING , or ATTENTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.


Classifications of Precautions

 DANGER	Improper use of this device poses a severe and likely risk of death or serious injury.
 WARNING	Improper use of this device may pose a risk of death or serious injury.
 ATTENTION	Improper use of this device may result in fire, injury, or property damage.




Symbols

	This symbol shows what you must not do.
	This symbol shows what you must do.
	This symbol shows what you must be careful with.

DANGER

 Prohibited	Do not use under any of the following conditions. Use of the device could cause an electronic medical device to malfunction.
	<ul style="list-style-type: none"> ◆ Persons using electronic medical devices. ◆ Pacemaker or other implanted electronic medical device. ◆ Artificial heart-lung machine or other electronic life-support medical device. ◆ Electrocardiograph or other electronic medical devices worn on the body.

WARNING

 Prohibited	<ul style="list-style-type: none"> • Do not use it with other beauty equipment at the same time. • Do not use it for purposes other than beauty, and do not use it in a damaged or malfunctioning state. • Do not look into or look directly at the RF LED lights. • Always move the device over the skin while treating. Do not use it on the same site continuously, Avoid damage to the skin or infection, inflammation, bleeding, pigmentation and other skin problems that will also lead to worsening of symptoms. • Do not use on body parts that can not be reached by the entire head surface. May cause a burn injury or a low-temperature burn injury. • Do not use the device on pets.
 Do not use near water.	<ul style="list-style-type: none"> • Do not use near water, in humid places, the bathroom or while bathing • Do not leave the device in a place with damp or water (on the washstand, etc.). Failure to comply may cause an electric shock, short circuit, or fire.
 Do not wash the device in water.	<ul style="list-style-type: none"> • Be sure that no water gets on the the unit or AC adapter cord. • Do not wash the device in water. Failure to comply may cause an electric shock, short circuit, or fire.



Do not use or touch with wet hands.

- Do not use or touch with wet hands.
Failure to comply may cause an electric shock, short circuit, or fire.



Prohibited

- Do not use the device near a flammable substance (such as alcohol, thinner, benzene, nail polish remover, or spray).
Use only with the supplied AC adapter, USB cable.
- Do not use the supplied AC adapter, USB cable for other devices.
Failure to comply may cause an electric shock, short circuit, or fire.



Disconnect the plug

- When unplugging the AC adapter from the unit or a power outlet, do so by gripping the plug.
- Do not unplug the AC adaptor by pulling.
- Do not use if the AC adapter is damaged or power socket is loose.
- In the following unlikely events that an abnormality arises, immediately unplug the AC adapter.
Disconnect the AC adaptor from the device and the plug before cleaning the device.
- Turn off the power and disconnect the AC adaptor from the power outlet when not in use.
Failure to comply may cause an electric shock, short circuit, or fire.



Careless or reckless use is prohibited

- Do not cut, damage, modify, excessively bend, pull, twist, or wrap the AC adapter.
Placing a heavy items on top of the AC adaptor or pinching them between other items, could damage them and cause fire or electric shock.



Must comply

- Use home-use power socket only.
- Make sure that the AC adapter is fully inserted into an electrical outlet.
Failure to comply may cause a fire, electric shock or accident.






Prohibited

- Do not connect the device to metal or electrical products (mobile phones, TVs, etc.)





User restrictions

- Keep away from children and do not allow children to use the device.
Failure to comply may cause an accident or injury.

 Disassembly prohibit	<ul style="list-style-type: none"> Do not disassemble, attempt to repair, or modify this device in any way. Failure to comply may cause a fire, electric shock or injury. Contact your local distributor for advice on repairing the device.
 Prohibited	<ul style="list-style-type: none"> Do not put in fire; Do not charge, use or leave in hot place.
 Must comply	<ul style="list-style-type: none"> Before using the device, please be sure to take off metal jewelry, glasses or contact lenses. Metal jewelry: earrings, necklaces, watches, etc. Always select Level 1 upon the first use until you become accustomed to use. Do not press the probe hardly on the skin and use continuously on the same body part for more than one minute, avoid skin problems. Close eyes before the treatments for the skin around the eyes starts. Stop using if you feel any physical discomfort during the treatments. Carefully check that you skin is showing no abnormality during a treatment. Do not use the device if a skin abnormality is perceived and consult your doctor for help.

WARNING

When using outside of China	
 Must comply	<ul style="list-style-type: none"> Be sure to use in an operating voltage range from 100 to 220 volts. Failure to comply may cause a fire, electric shock, or malfunction.
 Combustion warning	<ul style="list-style-type: none"> Be sure to check the plug configuration and the voltage of the power outlet before using. Failure to comply may cause the generation of heat, a fire, or a malfunction of the device. The plug configuration and voltage varies depending on the region and the type of accommodation. Always be sure to check for compatibility before using the device. For information about plug adapters, consult a travel agency or a consumer electronics retailer.

1.3 Contraindications:

The following persons should not utilize the device.

May cause an accident or other problems.





Prohibited

- Those who are under 18 years of age.
- Those who are in their pregnant, menstrual or breast feeding periods.
- Those who are intoxicated.
- Those with excessive myopia ($-6D$ or higher).
- Have severe concurrent conditions, such as cardiac disorders, have a pacemaker or internal defibrillator.
- High blood pressure, communicable diseases, malignant tumors, blood diseases.
- Those undergoing treatment by a doctor, those in a physical growth stage, those with acute medical conditions, tubercular conditions, in poor health.
- Currently have cancer or a history of cancer, especially skin cancer or premalignant moles; Have suspicious lesions within the treatment area.
- Have had any surgical procedure, laser resurfacing or deep chemical peeling within the treatment area within the past three months or are still undergoing a process to complete healing.
- Those who had an operation on their eye(s) or faces; who have a disease of the eye(s) or face; those who are currently undergoing treatment for a disease of the eye(s) or face.
- Have permanent fillers within the treatment area; Have metal implants in the treatment area, such as gold threads.
- If you have temporary dermal fillers within the treatment area, consult your physician before using the device. This does not apply to dental fillings or implants, which present no problem.
- Have had Botox administered in the treatment area within the past three months.
- Have any active skin condition within the treatment area, such as, but not limited to, sores, psoriasis, eczema or a rash; Have any disease that is stimulated by heat, such as, but not limited to, recurrent Herpes Simplex (within the treatment area).
- Those who take medication that induces drowsiness.

Do not use the device on various sites as indicated below and outside of specified sites on the face.

May cause an accident or other problems.

 Prohibited	<ul style="list-style-type: none">• Sites outside of specified sites• Area around the heart• Swollen moles, viral warts, on the Adam's apple, eyeballs, eyelids, lips, mucous membranes in the mouth and etc. (except when using Cool Mode).• Pathological pigmentation, moles, sites with viral warts.• Sites with muscle tingling and inflammation caused by sunburn.• Sites with pathological pigmentation caused by itching, fever, physical stimulation, etc.• Surgical sites, private parts, anus, area around Adam's apple, ear, hair removed site areas or sites have recently been shaved.• Body parts with cuts or that have undergone cosmetic surgery; inside the mouth; in the pubic region; on mucous membranes; acne; areas inflamed due to suppuration; areas with itchiness or inflammation.• Around the eye area (when not using the device with Eye Care or Cool Modes) RF LED Mode can be used over eyelids and around the eyes. When using the LED light on RF LED Mode on the eyelids and areas around the eyes, close your eyes and hold the device 5cm away from the eyes.• Do not look directly into the RF LED light.
	If you are not sure, consult your doctor before use.

2. WHAT IS THE TRI-ELF?

2.1 Definition:

TRI-ELF is the Home use beauty device combined with 3 replaceable probes in 9 modes, mainly focus on different skin problems with exact solutions.

T1 Probe-Naive ELF is also called Cleansing & Moisture probe, 5 modes.

T2 Probe-Mature ELF is also called Lifting & Tightening probe, 1 mode.

T3 Probe-Elegant ELF is also called Anti-aging probe, 3 modes.

2.2 Technology:

Use the TRI-ELF at home to help you solve the basic skin problems and reduce the appearance of mild to moderate facial wrinkles with different stage of skin problems in non-invasive treatment .

Science, Applied to Beauty. The TRI-ELF device is combined with 6 mainly technology. Such as radio frequency (RF), BIO, EMS, icon cleansing, LED and Cooling.

- Radio frequency (RF): stimulate collagen regeneration to achieve the results of anti-aging by creating heat energy in dermis.
- BIO: combination electroporation and electroosmosis technology by biological electricity to achieve the result of effective absorption for skin care products.
- EMS: electrical muscle stimulation is known as muscle electrical movement technology to achieve the result of lifting and tightening and outline shaping .
- Icon cleansing: absorbs dirt which cannot be removed with normal face washing with the power of ions.
- LED: LED Photo treatment for whitening & brightening, anti-acne & anti-inflammatory, water and oil balance etc.
- Cooling: soothing sensitive skin with cooling effect.

2.3 Story:

The TRI-ELF collections inspired by Triangel .

TRIANGEL, the guardian angel in Italian mythology, whose duty is to protect and guard the fate and help them find the most fate lover.

ELF is a type of humanoid supernatural being in Germanic mythology and folklore. In medieval Germanic-speaking cultures, elves generally seem to have been thought of as beings with magical powers and supernatural beauty.

BCD makes ELF means even more. E stands for EMS, L for LED therapy, and F for radio frequency.

Designed for our angels.

Naive angel in 25+, Mature angel in 35+, Elegant angel in 45+ Angel for you, elf for beauty, and

TRI-ELF for your exclusive beauty.

3.WHAT'S INSIDE THE TRI-ELF BOX ?

TRI-ELF Kit includes:

Device body*1

Electrode head*3

AC adapter*1

Make up Cotton pads*1

User manual*1



4.GETTING TO KNOW YOUR TRI-ELF DEVICE

4.1 General guidance before use:

Check the device condition before use

Always make sure the device in normal condition before operation.

(Check if blue indicator light can work normally or not, if there is any cracks and damage)

*Malfunctions may cause injury or accident.

Before Starting

- Remove all jewelry, earrings, glasses and contact lenses (if you wear them)
- Sit in front of a mirror during treatment so you can have better experience.
- Clean your face with facial cleanser and water and dry thoroughly.
- Make sure your hands are clean and dry.



How to apply gel

1. Dab a proper amount of skincare products on your fingertip from the surface of the preparation .
2. Dab proper spots of skincare products around the face.

△ Important

Be careful not to apply skincare products into your eyes.

3. Evenly spread a thin layer of the skincare products onto the face.

4. Clean your hands.

△ Reminder:

Please use the correct & recommended type of skincare product when you operate the different electrode head in different working mode.



How to connect the power supply and start the device


1. Install the probe
To positioning hole, Straight plug & lock, Rotation prohibited.
2. Connect the power supply cord to the device.
3. Plug the power supply into an electrical outlet.
4. Make sure the blue indicator light is on.

The device emits short beep along with indicator showing and blue light circle surrounded when ready to work.

5. Press the <M> button with <^> to turn the device on.

Remark:

- Adjust intensity level to ZERO with button (V) to turn the device off.
- The <LED> button only works under the T3 Probe-Elegant ELF (Anti-aging) in 3 colors with 5 Options.

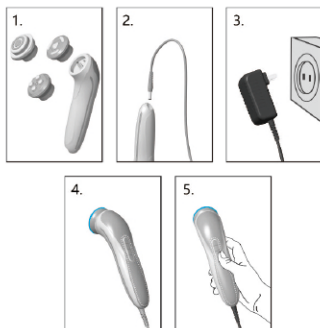

4.2 Product structure:

Handlepiece






Indicator panel

Level indicator:
adjustable 4 levels

Mode indicator:
optional 9 modes for
different skin care purposes



3 Replaceable Probes

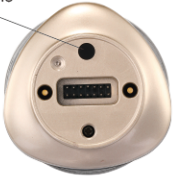	T1-Basic Skin Care Probe <5 Modes>	T2-Lifting & Tightening Probe <RF Mode>	T3-Anti-aging Plus Probe <3 Modes>
Naive ELF Probe	<p>LED lamps (Red, Orange, Blue)</p>  <p>Cotton stopper (Silver Ring) - Use this stopper to fix cotton to the electrode head.</p>	Mature ELF Probe	Elegant ELF Probe
		<p>Electrode head</p>  <p>Temperature sensor</p>	<p>LED lamps (Red, Blue, Green, Multi-Color)</p>  <p>Electrode head</p>
	 <p>TRI-ELF</p> <p>Intensity Indicator</p> <p>Mode Indicator</p>	 <p>Positioning hole</p>	

5. TREATMENT PROCEDURE

5.1 First time use

Pre-testing on skin

Before your first treatment, perform a test on a small area of skin (forearm, 2"x2"). This will ensure that your skin does not have a negative reaction (excessive redness and/or swelling) to the gel or treatment and that the device is operating correctly.

Apply a small amount of skin care products (check the types on the table of # 5.4) on the test area. Treat the area with the TRI-ELF. Follow the steps below:

1. Install the probe.
2. Make sure the blue indicator light is on.
3. Press the <M> button with <^> to turn the device on and adjust intensity level to ZERO with button (V) to turn the device off.
4. Place the device on the skin in right operations (check details on the table of # 5.5)
5. Check if your skin has a negative reaction (excessive redness and/or swelling).

Remark:

- If you experience any skin abnormality (rash, redness, itch, stinging or burning sensation) at the pre-testing site within 24 hours, immediately stop the pre-test and wash the pre-testing site and stop to use the device.
- If using your regular skin lotion, be sure that the skin lotion is water soluble.
- Do not use oil-based or scrub-based lotions, as this could cause the device to malfunction or skin problems.

5.2 Basic notice

- Be sure to move the unit on your skin while using it. Using the device continuously on the same area may cause burns or other injuries.
- The temperature of RF may feel different depending upon the condition of your skin and types of cosmetics used together.
- Do not use on eyelids or around the eyes (except when using the device on Eye Care or Cool Mode)
- Do not look directly into the RF LED light.
- Recommended ambient temperature is 15°C ~ 30°C
- Please take off the metal objects you are wearing (watch, necklaces, rings, earrings, etc.) before use. Otherwise, accidents and personal injuries may occur.
- In order to keep the skin moist, skin care products such as toner, gel and essence can be properly supplemented according to the skin condition during use.
Lack of moisture between probe and skin may lead skin problems because of decreased skin sensitivity.
- When the skin surface feels stinging, please lower the intensity level or change the use site. In addition, skin toner, gel, essence and other skin care products can be properly added to moisturize the skin, otherwise it may increase the skin tingling sensation.
- Do not use the device with or close to products containing alcohol, acid.
- Be used with or close to products containing alcohol, the applicator surface may be easily oxidized (discoloration, blackening, etc.), but it does not affect the penetration function.

5.3 Mode & level adjustment:

Select mode	T1 Probe-Naive ELF (Cleansing & Moisture)	5 Modes	Press (M) for Clean/Moisture/Eye Care/ EMS Up/Cool mode swift.
	T2 Probe-Mature ELF (Lifting & Tightening)	1 Mode (No adjustment)	Press (^) to work with RF mode, no need press (M) . * Mode cannot be selected.
	T3 Probe-Elegant ELF (Anti-aging)	3 Modes	Press (M) for MFEP/UP/EMS mode swift.
Adjust level	Clean/Moisture/EMS Up/ RF/MFEP/UP/EMS	4 Levels	Press (^) (V)for level adjustment.
	Eye Care/Cool	1 Level (No adjustment)	Press (^) to work. *Intensity level cannot be selected and no level display.
LED mode swift (Only work for T3)	MFEP/UP/EMS mode under T3 Probe-Elegant ELF	5 Options	Red/Blue/Green/Static 3-color/ Flashing 3-color

5.4 Treatment modes recommendations:

Probes	Product Description	Modes	Functions	Skincare solutions	Time	Intensity levels	Frequency (rough guide)
T1 Probe- Naive ELF (Cleansing & Moisture)	RF Ion Cleansing	Clean	Deep Cleansing	Cotton Toner	Approx. 3 minutes	4 levels	2 to 3 times a week
	RF LED(Red/Amber)	Moisture	Moisturizing	Essence Facial mask/1 piece	Approx. 3 minutes	4 levels	2 to 3 times a week
	RF MFEP	Eye care	Crow's feet removal	Eye essence	Approx. 3 minutes	4 levels	2 to 3 times a week
	EMS LED(Red)	EMS Up	Fine line removal	Essence Facial mask/1 piece	Approx. 2-3 minutes	N/A	2 to 3 times a week
	Cooling LED(Blue)	Cool	Soothing sensitive skin	Gel	Approx. 3 minutes	N/A	Daily use
T2 Probe- Mature ELF (Lifting & Tightening)	RF	RF	Lifting & Tightening	Essence Gel	Approx. 10 minutes	4 levels	2-3 times a week
T3 Probe- Elegant ELF (Anti-aging)	RF MFEP LED (3 colors with 5 Options)	MFEP	Lifting & firming plus	Essence Gel	Approx. 10 minutes	4 levels	2 times a week
	RF UP LED (3 colors with 5 Options)	UP	Effective absorption of nutrients	Essence Gel	Approx. 10 minutes	4 levels	2 times a week
	EMS LED (3 colors with 5 Options)	EMS	Deep wrinkle removal	Essence Gel	Approx. 10 minutes	4 levels	Daily use

Notes:

For the use frequency, please refer to the treatment method.

The use frequency and duration vary for each treatment head and each function mode.

Excessive care may cause adverse effects on the skin, please use according to the skin condition.

5.5 Use & operation method:

T1 Probe-Naive ELF
Cleansing & Moisture



Use Method


Attaching Cotton to the Head

1. Remove the cotton stopper from the electrode head.
2. Place the cotton on top of the electrode head and press the cotton stopper down over the cotton.
3. If the cotton stopper does not fit properly on the electrode head, use a thinner cotton.
4. If the cotton is too thick, the thermal action of RF treatment may be adversely affected and you will feel less heat.
5. Apply your regular skin toner onto the cotton.

* Make sure that the toner coats the entire surface of the treatment head.

* This device is not water-proof so be sure to thoroughly wipe off any toner left on the device with a dry cloth after each use.

Do not wash with water or other liquids.

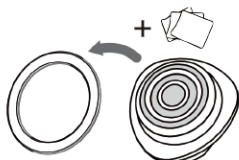

Operation Method

<Clean Mode & Moisture Mode>

Move the probe more quickly when the temperature is overheated shortly. In addition, the heat will be reduced with the use of make up cotton. Or You can adjust to a comfortable intensity level using the Level button. Always move the device over the skin while treating. Do not use on body parts that cannot be reached by the entire probe surface.



⚠ CAUTION

No treatment area:

Eyeballs, eyelids or eye area,
Adam's apple.

Recommended sliding speed:
2cm/s

Left face: Approx.
1.5 minutes

Right face: Approx.
1.5 minutes

<EMS Up Mode>

When the skin surface feels stinging, please lower the intensity level or change the use site. In addition, the gel or essence must be properly added to moisturize the skin, otherwise it may increase the skin tingling sensation.



⚠ CAUTION

No treatment area:

Eyeballs, eyelids or eye area,
Adam's apple.

Circle movement as indicated
by the arrows.

Left face: Approx.
1.5 minutes

Right face: Approx.
1.5 minutes

<Eye Care Mode>

Apply gel or eye essence and operate the device from the corner of your eye to the end of your eye, and repeat the steps as press-release-move.

*The electrode head cannot be operated at a fixed position more than 3 seconds.

*Using slight pressure ensure that energy tips firmly touch the skin at all times during treatment.



CAUTION

No treatment area:
Eyeballs, upper eyelids or eye area, Adam's apple.
*lower eyelids, the corner of eye can be operated within 3 seconds.

Left face: Approx.
1 minutes

Right face: Approx.
1 minutes

<Cool Mode>

Apply gel and operate the device with steps as press-release-move. (repeat)

CAUTION

No treatment area:
Eyeballs, upper eyelids or eye area, Adam's apple.
*lower eyelids, the corner of eye can be operated.

Left face: Approx.
1.5 minutes

Right face: Approx.
1.5 minutes

T2 Probe-Mature ELF
(Lifting & Tightening)



Use Method

Apply an even layer of gel or essence to the whole treatment area (face or neck). During the operation process, appropriate gel or essence can be added according to the skin condition. Always select Level 1 upon the first use until you are accustomed to use.

Operation Method

<RF Mode>

1. Identify the parts of your face that you wish to treat.
Divide your face into three treatment areas.



2.TREATMENT SAFETY GUIDELINES

Using slight pressure ensure that all four energy tips firmly touch the skin at all times during treatment.

It is important to keep the device moving at all times during treatment.

Leaving the device on the same spot without movement may result in skin overheating and possible excessive skin redness.

For most effective treatment, move the device slowly in circular motions.

circular movements reference



3.about 3-5 minutes Slow and circular movements for each zone, Approx.10 minutes for whole treatment.


For most effective treatment, using lifting movement from jaw line to cheeks.

Security Mechanism

The security mechanism can protect your skin from accidents or a injury caused by over-heating. The device will emit beep sounds and automatic power off at a over-heating temperature. The device will restart and stop beeping when the temperature decreases to a safe temperature, or the device has been removed to a substance with low temperature.



Never treat directly on eyes / eyelids, lips, nose, ears or thyroid.

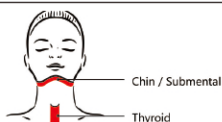

■ Recommended treatment areas

■ Areas prohibited for treatment

Do not use if you have moderate or severe wrinkles underneath the chin

T3 Probe-Elegant ELF
Anti-aging



Use Method

Apply an even layer of gel or essence to the whole treatment area (face or neck). During the operation process, appropriate gel or essence can be added according to the skin condition. Always select Level 1 upon the first use until you are accustomed to use.

Operation Method

<MFEP Mode & Up Mode>



Lifting operation with slow motion, 2-3 times for each step.

*The electrode head cannot be operated at a fixed position.

√During the treatment you must keep the device constantly moving.

√If your skin begins to feel too hot, expand the treatment area or move slightly faster.

*During the operation process, appropriate gel or essence can be added according to the skin condition.

<Ems Mode>



Apply the gel or eye care essence on the skin, then apply the device as if you are "stamping" onto your skin, Repeat the steps as press-release-move. 2-3 times for each step.

*Do not press the electrode head against the skin for longer than 3 seconds.

*During the operation process, appropriate gel or essence can be added according to the skin condition.

The treatment time for full face skin care is about 10 minutes. priority to key sites. The treatment time varies for different users.

5.6 LED:

T1 Probe-Naive ELF

Cleansing & Moisture



Moisture Mode
LED (Red/Amber)

Repair and improve skin barrier



EMS UP Mode
LED (Red)






Accelerates micro-circulation for shiny and elastic skin



Cooling Mode
LED (Blue)

Oil & Water balance, anti-acne and anti-inflammatory

T3 Probe-Elegant ELF
(Anti-aging)

	Static 3-color	complex skin troubles or for anti-aging skin care.
	LED (Red)	tired skin due to aging; whitening & brightening
	LED (Blue)	acne and other skin problems. anti-acne & anti-inflammatory
	LED (Green)	dull & sensitive skin water & oil balance, soothing redness.
	Flashing 3-color	enhance the essence penetration and absorption in the last procedure of the treatment.

6. CLEANING AND STORAGE

If you will not be using the device for an extended period, clean the device before storing.

How to clean and store the product

- When the device is not in use, turn off the main body and disconnect the AC adapter from the main body and wall outlet.
- To clean the electrode head, wipe it with a dampened and well-wrung cloth to remove dirt. Wash the cotton stopper with water.
 - Remove the cotton stopper from the main unit, and wash it with water then wipe dry.
 - Refrain from using an alkali cleaning agent, alcohol, thinner, benzene, nail polish remover, or other such substance.
 - Do not leave the unit with lotion or gel on it.
- Wipe off any dirt from areas other than the electrode head and cotton stopper with a dry cloth.
- Avoid storing the device near a stove or heater or in a bathroom or other high-temperature/high-humidity location, location where there is exposure to direct sunlight, or unstable location.
- Keep the device out of the reach of children or pets.
- Do not clean the device with any wet wipes. (Usually the wet wipes contain alcohol ingredient.)

7. TECHNICAL SPECIFICATIONS

Specification	
Description	TRI-ELF/3 in 1 Home Use Beauty Device
Power	20W
Input Voltage	AC110/220V, 50/60HZ
Color	Gold
LED Wavelength	Blue: 415nm (± 10 nm) Green: 530nm (± 10 nm) Red: 620nm (± 10 nm) Yellow: 590nm (± 10 nm)
Frequency	T2-1MHz; T3-94K
Package size	26*22*9 CM / VWK: 1.03 kg
G.W.	0.85kg

8. WARRANTY TERMS

8.1 Warranty terms

This product has been subjected to strict quality control performed by the manufacturer. For this reason, the distributor in your country guarantees proper function of your product for 12 months from date of actual purchase. If you use the device properly, follow all instructions outlined in the user manual, you can enjoy the product to your ultimate satisfaction. This policy warrants the support that the distributor provides for the device, should it happen that failures attributable to manufacturing defects be found. In such case, you shall be entitled to make full use of the benefits of this warranty. This implies the reason and replacement of parts, free of charge (regarding both service and parts) at any Authorized Service Center. Within its validity term, you are entitled to make use of the warranty if necessary to have a fully enjoyable product, provided the above-indicated conditions are fully met. This warranty is valid for the device marketed by the distributor in your country.

The warranty is not valid when:

1. Amendments in the warranty policy, sales receipt or invoice, or the absence of any of these original documents.
2. Damage or deterioration of the device due to improper use or if the device has been modified or repaired by an unauthorized third party.
3. Damage of the product due to defects incurred during transportation, such as bumps, bangs or improper handling.
4. Damage due to humidity or any kind of liquids, foreign bodies inside the device or the use of unapproved, unadvised accessories with the device contradictory to the user manual.
5. When the device is not used in accordance with the respective user manual or with its supplied accessories.
6. Partial or total damage or loss of the device, due to natural disasters (such as earthquakes, floods, fire, lightning, etc.), force majeure or accident.
7. Damage to the device after being subjected to high dust environments, humidity or application of excessive voltage from the power inlet.
8. To cosmetic damage, including but not limited to scratches and dents.

8.2 Warranty card

The right of interpretation belongs to the brand or manufacture.

Commodity	3 in 1 Home Use Device
Model No.	TRI-ELF
Customer Name	
Address	
Email	
Warranty Start Date	
Warranty End Date	

We hereby guarantee and warranty all products owned for the period of time.

Make Beauty Easier