

装
订
线



BB Body User Manual

TABLE OF CONTENTS

1.PRODUCT OVERVIEW.....	02
2.INSTRUMENT PARAMETERS.....	02
3.PRODUCT DIAGRAM.....	03
4. OPERATION INSTRUCTIONS.....	04
5.PRECAUTIONS.....	07
6. PRECAUTIONS.....	07
7. CONTRAINDICATIONS.....	08
8.MAINTENANCE AND CARE OF THE DEVICE.....	09
9. WARRANTY CARD.....	10

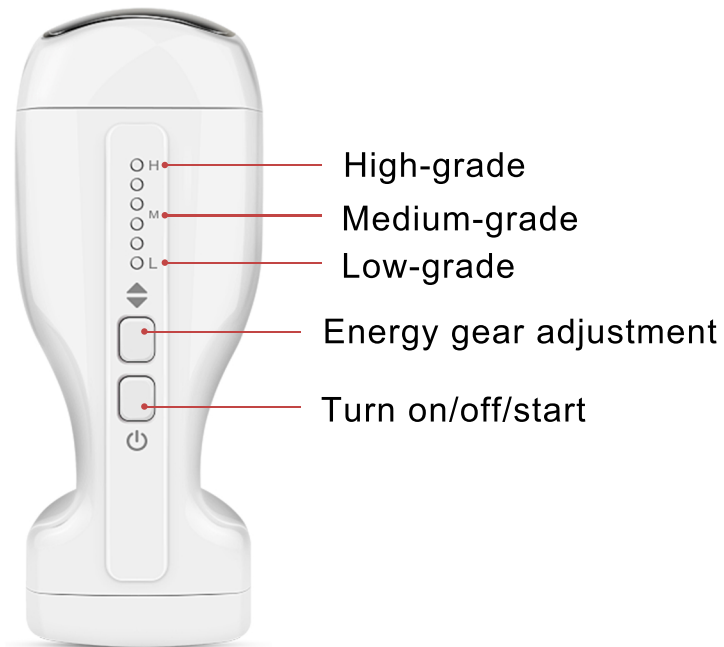
1.PRODUCT OVERVIEW

Welcome to use BB Body, a versatile beauty instrument designed for both home and salon use. Before using the machine, please carefully read this manual and strictly follow the instructions to operate this beauty device.





2.INSTRUMENT PARAMETERS

Product name	Bb Body	Power supply voltage	100-240V
Product model	H-038	Product Dimensions	7.6*18.5*7.8cm
Working current	3A	Working voltage	DC24v
Power	72W	Treatment Heads	1

3.PRODUCT DIAGRAM




4.OPERATION INSTRUCTIONS

- 1.Before using the BB Body, ensure that the treatment area is clean and dry.
- 2.Insert the power cord into the power outlet and press the power button for 3 seconds to turn on the device . The BB Body will be powered on.
- 3.Once the device is powered on, lightly press the energy adjustment button  to select the desired energy level. Start with the lowest level (L)  and gradually increase the energy level as needed. To start the energy output , press the start button.

Note: To ensure the longevity of the device, it is equipped with an automatic shutdown feature after 30 minutes of continuous operation. If you need to continue using the device, simply follow the standard startup procedure.

Specific Operation:

Apply a suitable medium (gel, aloe gel, essential oil, fat-burning milk, slimming radiofrequency cream, or other products) to the treatment area.Place the handle on the treatment area and lightly press the handle's start button . The handle will display a green or red light (green light indicates low energy, and the light will gradually turn red as the energy level increases). Please be mindful of the energy level and immediately lower it if the sensation becomes uncomfortable. The device will start energy output.

Note: During usage, you may experience a warm sensation on the skin and a buzzing sound in the ears. This is normal. The device uses strong sound waves at 40,000Hz to penetrate the subcutaneous fat cells, causing resonance and producing powerful shock waves and free radicals within the fat cells. This effectively breaks down fat cells into fragments and free fatty acids.

Operation Time Intervals:

For arms and thighs: Treat every 2 days for the first 10 sessions, then treat every 3 days after the initial 10 sessions.

For the abdomen: Treat every 3 days for the first 10 sessions, then treat every 5 days after the initial 10 sessions.

Note: During the treatment, remember to perform slow, continuous movements and avoid stationary treatment for extended periods, as it may cause skin damage.

During the operation, ensure the treatment area remains moist. When using the device's metal probe, rotate and maintain constant contact with the skin by alternating the four corners of the probe. If only two or three corners are used, the untreated corner may experience higher temperatures, which could potentially damage the skin. Always ensure all four corners are rotating during the operation. Please adhere to these guidelines for effective and safe operation.

Thigh Operation:

Treat each thigh for 20-30 minutes. During the operation, move or make circular motions from the outer to inner thigh. Then, slide and lift from above the knee towards the base of the thigh.

Abdomen Operation:

For the abdomen, it is recommended to treat for 30-45 minutes. For individuals with more body fat, the treatment time can be extended to 45 minutes, while those with less fat can treat for 30 minutes. During the operation, move in circular motions or make linear movements from the side waist towards the navel. Then, perform clockwise circular motions near the navel area.

Arm Operation:

For each arm, the recommended treatment time is 15-20 minutes. Perform movements from the outer to inner side of the arm and lift from bottom to top. Focus on areas with excess fat, such as the "bye-bye" arm area.

Please follow these guidelines for effective and safe treatment.

5. PRECAUTIONS

- 1.Avoid using pasty substances; prefer essential oils as a medium.
- 2.The device will automatically shut down after a period of use (approximately 30 minutes). To continue using it, restart the device.
- 3.After treatment, in the final step, you can use massage techniques to alternate wiping and draining towards the nearby lymph nodes to promote better toxin metabolism.
- 4.For loose fat, the treatment time can be extended by approximately 5 minutes.
- 5.It is recommended for customers to drink a glass of warm water before the treatment, as it helps to accelerate skin blood circulation and eliminate toxins from the body.

6.PRECAUTIONS

- 1.It is recommended to drink at least 2000cc of water daily (equivalent to eight glasses).
- 2.Reduce consumption of spicy foods, smoking, alcohol, and late nights. Increase the intake of fruits and vegetables and opt for light meals.
- 3.Avoid high-temperature environments like saunas for one week after treatment.

7.CONTRAINDICATIONS

- 1.Avoid treating areas with skin wounds or scars.
- 2.Pregnant women and women during menstruation should not use this device.
- 3.The device should not be used for individuals with internal metal implants, such as pacemakers.
- 4.People with extremely weak physical conditions should only use the device after improving their health.
- 5.Women who have inserted an intrauterine device should avoid treating the abdomen within three months.
- 6.Women who have given birth within three months or undergone a cesarean section within six months should not treat the abdomen.
- 7.Individuals with skin allergies should not use any beauty device during an allergic period.
- 8.Individuals experiencing allergies are advised to wait for a week after the allergic reaction subsides before using the device.
- 9.Avoid treating areas where implants or metal are present.
- 10.Avoid treating areas with wounds or large abscesses.
- 11.Patients with serious illnesses, heart disease, diabetes, mental illnesses, etc., should use the device with caution.
- 12.Patients with severe skin diseases or infectious diseases should use the device with caution.

8. MAINTENANCE AND CARE OF THE DEVICE

1. After each use, clean the device with water and store it properly.
2. Before using the device, ensure that it is plugged into an earthed power socket.
3. Ensure the device's voltage is suitable. In case of unstable local power supply voltage, use a stabilizing power supply for compatibility.
4. To maintain the treatment effect and ensure the normal service life of the device, use only the specified accessories provided or recommended by the original manufacturer.
5. Keep the device away from damp areas or water sources, and avoid direct exposure to sunlight.
6. Avoid placing the device near strong heat sources, as it may affect its lifespan and proper functioning.
7. Do not allow water, oil, or other substances to enter the internal components of the device. Avoid hitting or dropping the device to prevent damage.
8. The device comes with a one-year warranty. Within this period, if the machine experiences non-human-induced damage, the manufacturer will provide free repairs and replacement of parts (round-trip shipping costs to be borne by the customer, while the factory will cover the replacement of parts and repair costs).
9. After the warranty period ends, the manufacturer will provide lifelong maintenance for the device (with the manufacturer charging only the cost of replacement parts and not charging any repair fees).

9. WARRANTY CARD

Warranty Card	
Product name:	Model:
Purchasing date:	Store name:
Address:	Tel:
Customer name:	Wayyanty period:
Failure cause:	
Repair content:	
Warranty terms:	
1. If the product has quality problems caused by non-human damage, it will not be replaced if there is no quality problem.	
2. If it is a software failure, the machine can be repaired within three days without replacement.	
3. If the product has hardware failure caused by non-human factors, the company will provide warranty within one year from the date of purchase of the product.	
4. The machine does not turn on, no warranty.	
Remark:	