

Professional Picosecond Device

Instruction



LS-008



Content

Importante attention matters.....	01
Safety attention matters.....	01
Button introduction.....	03
Operate way.....	04
How to use repair cream.....	06
Using steps.....	06
Operation methods of different affected areas...	08
Skin care.....	10
Operate suggestions.....	11
Specifications.....	12
Simple fault exclusion table.....	13
Warranty Card.....	14

Thanks very much for choosing our pico second laser pen. Before using the product, pls read and understand all the instructions and precautions.

Importante attention matters! ⚠

When using this product for the first time, suggest to find something with color to practice(for example black spots on fruit or paper,etc.) After being proficient in the operation specifications, you can operate on your body. Recommended to start from the arm and finally operate on the face or other parts.

For the first time, it must be used from the lowest gear, and the gear can be gradually increased according to your actual conditions. Each person's skin sensitivity is different. Please do not rashly increase the gear.

Safety attention matters! ⚠

Taboo

- ① Child and elderly;
- ② Light -sensitive skin and people who are using photosensitive drugs;
- ③ People who input heart pacemaker and other medical electronic machines


Attention:

- ① Women during pregnancy, or prepared to pregnant, breastfeeding women are not recommended to use;
- ② It's forbidden to aim the laser light outlet window at the fragile parts such as eyes,ect.,
- ③ Take the eyes protector glasses when using it;
- ④ Even if you wear the eye protector glasses,do not look at the laser light outlet window directly;
- ⑤ It is forbidden to operate devices in a flammable and explosive environment; keep the room well ventilated;
- ⑥ Please operate at the place where children are far away, and place it in a place where children cannot contact after use;

Prohibit:

- ① Do not transform, disassemble or repair this product by yourself.
- ② Do not use this product in a flammable and explosive environment.

Button Introduction

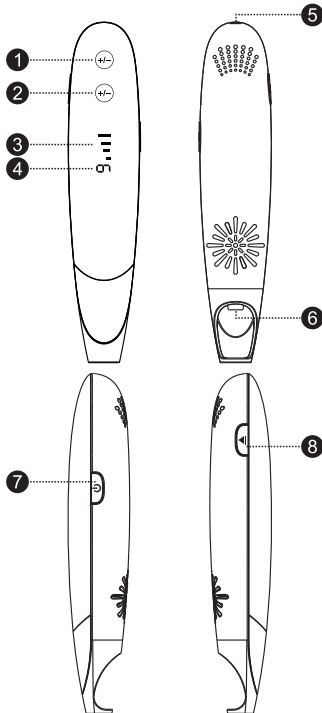

- ① Flash intensity adjustment button
- ② Flash speed adjustment button
- ③ Light flash speed display: have 4 gears, 1 gear is the the slowest flash speed , 4 gear is the fastest flash speed , 1 gear is default when turned on.
- ④ Light flash intensity display: have 9 gears, 1 gear is the the minimum intensity , 9 gear is the largest strength, 1 gear is default when turned on.
- ⑤ Power interface
- ⑥ Laser exit
- ⑦ Power on/ off/ launch button: press it for 2s will turn on/off, press the "location" button, don't move, and then press "launch" button, it will shot laser.
- ⑧ Location button: press it when power on, it will shot blue guide lighting, accurately locate the affected area.

Operate way

1. Start and test the machine

Plug in the power supply, the machine makes a "di" sound, and the light intensity adjustment button lights up

2. Gears selection

Remove freckles, moles, tattoos, etc., the gears selected is flexibly adjusted according to its own situation and skin area. Those who are skilled can use high -gears, fast speed, and convenient operation. When using it for the first time, you must adjust the intensity and speed of the laser to 1 gear, and use the minimum intensity. You can gradually increase the level according to your own situation. It can choose flexibly. Please do not rashly increase the gear!

3. Aiming at the affected area

In actual operation, the best distance is skin touched by the alloy bracket. The center of the metal ring on the head of the alloy bracket is the laser point. Place the effected skin in the center of the metal ring, and press the "location button" for aiming.

4. Wear glasses

Press the "location button" without remove, then press the "Launch button". After Didi hit a few times, you can release the button and open your eyes. Repeat operation in this cycle, you can treat the affected area until it in best standard. (See: The operation methods of different affected areas.)

5. Do not aim at the eyes

Laser will act to dark places. Do not let minor touch it. Do not aim at your eyes.

How to use repair cream

Repair paste

After using, keep dry and hygienic, let the affected area be naturally scab. After the skin is completely scab on the second or third day, you can use repair cream twice a day. You can also use repair paste, scald cream, moist cream and other professional ointment, use according to its instructions.

Using steps


1. Wear protected sunglass.




2. Plug in the power supply, the machine makes a "di" sound, and the "flash intensity adjustment button" lights up.




3. Press the "Power on" 2s to turn on

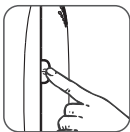


4. Press the "Flash speed adjustment button" in short time to choose flash speed;

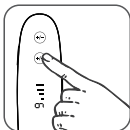

5. Press the "Flash intensity adjustment button" in short time to choose flash intensity;

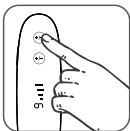

6. Press the "location button" laser outlet will emit guidance light, located to the skin therapy area, don't move, press the "Launch button" to shot laser.

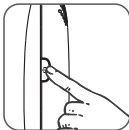

Attention: Only press the "location" and "launch" buttons at the same time can emit laser.

7. After using, press "Power off" button will shut down the device, will shut down automatically if no operation in five minutes.

8. Please clean it after using, and then safe keeping

Operation methods of various affected areas

1. Remove speckle

Freckles, genetic spots, pigment spots, elderly spots, sliced spots, etc., can be all removed at one time. Aim the product at the location of the spots and shot the spots until it become red or yellow, and the spots will fall off with the crusting in the following days. In addition, the skin spots on face are not recommended to operate. If there not new spots growing within two years. It is in stable period and that can be used.

2. Washing eyebrows and tattoos

Eyebrow/tattoo requires a course of treatment. Generally need 3 courses, and one month is a course treatment. In

addition, the eyebrows and tattoos do not need to be wiped off when using. Every time you used on affected area, the epidermis become gray-yellow or the epidermis has a layer of crust. Using once a month, every time have a thin crust, the eyebrows and tattoos will perfectly washed after several times. If the affected area is too large, it should be processed multiple times. Do not complete it at one time, and accidentally suffer from bacterial infection.

3. Remove mole

Acne marks, pigment layers, moles can be wiped when using device. Use a cotton swab to dip physiological saline or rub it with alcohol cotton. You can use it until without melanin. In addition, it is necessary to pay attention to the depth of moles. Some moles grow on the surface of the skin, while the root of moles is deeply tissue.

4. Acne marks, pigment layers

The method of acne, moles and freckle removal are different. We need to learn to distinguish whether it is acne marks or moles. Acne marks are grow acne for a long-term or marks after acne scar falling off for months, this can disappear themselves, don't need to remove it. .

The marks keep for 1~2 years are moles, these can be removed.

5. Other parts

Because it is a picosecond laser, it will only act to the area where have melanin. For the affected parts of the light colors (such as: light color spots or marks, light color tattoos, claw, wart, etc. without melanin), you can use eyebrow pencils to blackened and then use, to make it have the best efficacy

Remind:

Each person's skin sensitivity is different. Some people may have blisters when use. Don't worry about it, that's doesn't matter. The main substances of the blisters are electrolytes, glucose, fiber protein, etc. The foam skin is a good protective film, it's had better to make it absorb it by itself. For big blisters, you can pierce it after with acne needles disinfection, and then use professional ointment. Do not touch the infection with your hands.

Skin care

1. Wash face and take a bath

Moles, small amounts of spots, do not need to worry about it. If it's small number, you can use a waterproof band aid to paste it temporarily, and remove it after washing. If it's much spots, and you cannot paste it with a waterproof band aid, you can use a warm towel to gently rub the face in the following two days after using,. On the third day, you can quickly wash your face with water, then dry with a dry towel gently. When the scab fall off, you can wash your face with facial cleanser.

2. How to take care

Can't drink before the scab falls off, and other can eat casually;

If the area is large, and causes infection accidentally, you can buy Iodine to wipe, and then apply the professional medicines such as redlyticin ointment, Scald Palace, etc., or follow the doctor's advice. In addition, let scab go off naturally, don't pull it with your hands!

3. About skin become red after using

There are generally red marks after crusting falling. Don't worry. This is common for skin growth. It is normal, that's the color of the new skin! Generally speaking, it will disappear by itself. If you want it to disappear faster, you can massage it. Take a small amount of vitamin E or aloe vera gel in the red seal. Use your fingers to rub it with your fingers.

Operation suggestions

When used for the first time, you may have various problems and concerns. This is normal, mainly because it has not been operated. It is recommended to practice with a colorful thing, such as black spot on fruits. After practicing, start with your body arm, and finally operate on the face, step by step, so that you have confidence in yourself! If you don't understand, you can contact the customer service, and the customer service will guide you.

Specifications

Model	LS-008
Size	210x45x38mm
Rated power	< 8W
Adaptor	Output: AC100-240V, 50/60Hz Input: DC5V --- 2A
Work Tem.	Temperature: 5 ~ 40°C Relative humidity: 20 ~ 90%, no dew
Storage Tem.	Temperature: -20 ~ 55°C Relative humidity: 5 ~ 95%, no dew

Simple fault exclude

Performance	Reason	Deal way
No power	The power cord is not connected successfully	Connect to the power cord
Cannot turn on	1. Didn't press the "Power" button with force 2. Didn't press the hold the "Power" button exceed 2s.	Press the "Power" button with a little exceed 2s.
No laser output	Only press the "Launch" button will not shot laser	Only when press the "launch" laser and "location" laser at the same time it will shot laser.
Shut down automatically	The device will shut down automatically if no operation in 5 minutes.	Press the "power" button for 2s to turn on



Warranty Card

Dear product user:

First of all, thank you very much for choosing our products, so that we have the opportunity to provide you with the best quality service. In order to make our services more satisfied, please read this instructions carefully and keep this warranty properly after buying Card. (Effective only within the People's Republic of China).

Warranty

During the warranty period, the protection of the products (for free maintenance) is started from the date of purchase. The products during the warranty period are only maintained at the brand identification location; when maintained, please show this warranty card, purchase time and voucher. From the date of purchase, during the following period, it is repaired under normal use; the whole machine warranty within 12 months (except commercial products) belong to one of the following situations, even during the warranty period, will not repair free of charge, but can repair in charge.

- Product damage caused by violating the instructions of the product use;

- Product damage caused by components without legally certified (such as chargers, etc.);
- Product damage caused by exceed the allowable use environment (such as high temperature, corrosion, etc.);
- Product damage caused by non -normal external forces (such as falling, squeezing, bumps, impacts, etc.);
- Products damage caused by natural disasters or other force majeure (such as fire, flood, earthquakes, etc.);
- Products that counterfeit and illegal smuggling;
- Cannot provide valid purchase vouchers, invoice, warranty cards, etc.;
- Damage to products or parts caused by other outer factors or improper use.

In order to provide better service, pls fill in the below info.:

Name: _____

Wechat: _____

Buy ID: _____

Contact: _____

Warranty Reason: _____

装订方式：骑马订装

尺寸：54×120mm

材质：封面157g铜版纸

内页105g铜版纸

不用过膜

封二封三空白，单色印刷。