


MOOLWEEL



MOOLWEEL

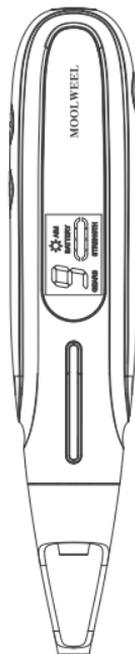

Picosecond laser beauty meter

**Instrument**

LS-006

## Content

Importante attention matters.....	01
Safety attention matters.....	01
Button introduction.....	03
Operate way.....	04
How to use soothing cream and repair cream....	06
Using steps.....	07
Operation methods of different affected areas...	08
Skin care.....	10
Operate suggestions.....	11
Specifications.....	12
Simple fault exclusion table.....	13
Warranty Card.....	14

Thanks very much for choosing the MOOLWEEL brand. Before using the product, pls read and understand all the instructions and precautions.

### Importante attention matters!

After being proficient in the operation specifications, you can operate on your body. Recommended to start from the arm and finally operate on the face or other parts.

For the first time, it must be used from the lowest gear, and the gear can be gradually increased according to your actual conditions. Each person's skin sensitivity is different. Please do not rashly increase the gear.

### Safety attention matters!

#### Taboo

- ① Child and elderly;
- ② Light -sensitive skin and people who are using photosensitive drugs;
- ③ People who input heart pacemaker and other medical electronic machines


#### Attention:

- ① Women during pregnancy, or prepared to pregnant, breastfeeding women are not recommended to use;
- ② It is forbidden to put the bald head at the fragile parts such as the eyes, ect.,
- ③ Take the blinkers when using it;
- ④ Even if you wear a blinkers,do not look directly at the lamp head light or reflection light;
- ⑤ It is forbidden to operate instruments in a flammable and explosive environment; keep the air in the room smooth;
- ⑥ Please operate at the place where children are far away, and place it in a place where children cannot contact after use;

#### Prohibit:

- ① Do not transform, disassemble or repair his product by yourself.
- ② Do not use this product in a flammable and explosive environment.

## Button Introduction

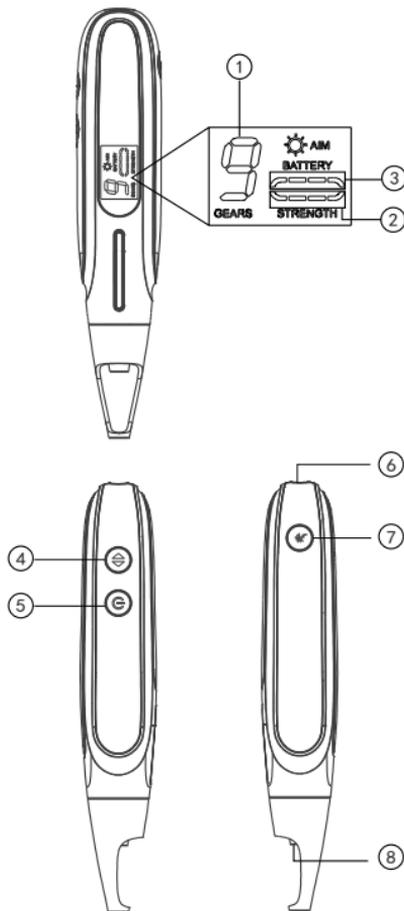

- ① Laser speed display: 9 gear in total, 1 gear is the fastest, 9 gear is the slowest, 1 gear is default when turned on.
- ② Laser intensity shows: 4 gear in total, 1 gear is the the minimum intensity , 4 gear is the largest strength,1 gear is default when turned on.
- ③ Battery electricity indicator
- ④ Laser strength adjustment button
- ⑤ Power/shut down/flash button/laser speed adjustment button: press for 2 seconds will switch between power on/shut down. Press the "Positioning button" without remove, and then press the "Launch Press"Button will shot laser. The short pressing of the boot is "Laser speed adjustment button" when power on.
- ⑥ Power interface
- ⑦ Positioning button: Press the "Positioning button" laser exit will send a guidance light to locate the skin.
- ⑧ Laser exit

## Operate way

### 1. Turn on the instrument

Press the "power switch" for 2s will turn on.

## 2. Gears selection

Remove freckles, moles, tattoos, etc., the gears selected is flexibly adjusted according to its own situation and skin area.

Those who are skilled can use high -gears, fast speed, and convenient operation.

When using it for the first time, you must adjust the intensity and speed of the laser to 1 gear, and use the minimum intensity. You can gradually increase the level according to your own situation. It can choose flexibly. Please do not rashly increase the gear!

## 3. Aiming at the affected area

In actual operation, the best distance is skin touched by the alloy bracket. The center of the metal ring on the head of the alloy bracket is the laser point. Place the effected skin in the center of the metal ring, and press the "Positioning button" for aiming.

## 4. Wear glasses

Press the "Positioning button" without remove, then press the "Launch button". After Didi hit a few times, you can release the button and open your eyes.

Repeat this in this cycle, you can hit the affected area in best stander. (See: The operation methods of various affected areas.)

## 5. Do not aim at the eyes

Laser will act to dark places. Do not let adults touched.

Do not aim at your eyes.

## How to use soothing cream and repair cream

### 1. Soothing cream

Clean the affected area first, and the area of soothing cream can be covered through affected area, and the thick around 1 coin. Then put a layer of plastic wrap on the outside (it will be dry without plastic wrap, that makes the effect is not good) about 20 minutes to 40 minutes (afraid of pain for 40 minutes.), Then scrub the soothing paste to start therapy.

(Tip: Soothing cream is strictly prohibited to drip into the eyes and in the mouth: otherwise pls clean your skin with cleaning water immediately and then seek medical treatment; tighten the refrigerator with endless use

### 2. Repair paste

After using, keep dry and hygienic, let the affected area be naturally knotted. After the skin is completely infatuated on the second or third day, you can use repair green twice a day. use repair paste, scald cream, Wet cream and other professional ointment, use its instructions.

## Charger

1. Put the plug of the charging adapter on the socket.
2. Connect the charging socket on the instrument. The electricity indicator symbol circulation light that it starts charging.



## Using steps

1. Wear protected sunglasses



2. Press the "Power" 2s to turn on



3. Press "strength button", select the strength



4. Press laser speed button, select the shot speed



5. Press the "Positioning button" laser outlet will emit guidance light, position it to the skin therapy site, don't let go, press the "Launch button" output laser.



**Attention:** Only press the "positioning" and "transmit" buttons at the same time can emit laser.

6. After using, press "Power" button will shut down the instrument, will shut down automatically if no operation in five minutes.

7. Please clean it after using, and then safekeeping

## Operation methods of various affected areas

### 1. Remove speckle

Freckles, genetic spots, pigment spots, elderly spots, sliced spots, etc., are all removed at one time. Aim the product at the location of the spots and shot the spots to red or yellow, and the spots will fall off with the crusting in the following days. In addition, the skin spots on face are not recommended to operate. If there are not new spots forming within two years. It is in a stable period and that can be used.

### 2. Washing eyebrows and tattoos

eyebrow/tattoo requires a course of treatment. Generally need 3 courses, and one month is a course treatment. In

addition, the eyebrows and tattoos do not need to be wiped off when washing. Every time you used on affected area, the epidermis is grayyellow or the epidermis has a layer of crust. Washing once a month, every time have a thin crust, the eyebrows and tattoos are perfectly washed in batches. If the area is too large, it should be processed multiple times. Do not complete it at one time, and accidentally suffer from bacterial infection.

### 3. Remove mole

Acne marks, pigment layers

Moles can be wiped when using device. Use a cotton swab to dip physiological saline or rub it with alcohol cotton. You can use it until without melanin. In addition, it is necessary to pay attention to the depth of Moles. Some moles grow on the surface of the skin, while the root of moles is deeply tissue.

### 4. Acne marks, pigment layers

The method of acne, moles and freckle removal are different. We need to learn to distinguish whether it is acne marks, long-term gangrenes or sinks mark within a few months, this can disappear itself, no need to remove it. You need to use for these marks more than one year.

### 5. Other parts

Because it is a picosecond laser, it will only reflect to the area where have melanin. For the affected parts of the light colors (such as light spots or marks, color light tattoos, meat grains, chicken eyes, rackets, etc. without melanin), you can use eyebrow pencils to blackened and then use, to make it have the best efficacy.

### Remind:

Each person's skin sensitivity is different. Some people may have blisters when use. Don't worry about it, that's doesn't matter. The main substances of the blisters are electrolytes, glucose, fiber protein, etc. The foam skin is a good protective film, it's had better to make it absorb it by itself. For big blisters, you can pierce it after with acne needles disinfection, and then use professional ointment. Do not touch the infection with your hands.

## Skin care

### 1.Wash

Moles, small amounts of spots, do not need to worry about it. If it's small number, you can use it with waterproof and paste it temporarily. Remove it after washing. If it's much spots, and you cannot stick to the waterproof and bandage. In the follow two days after using, you can use a warm towel to gently rub the face. On the third day, you can quickly wash your face with water, then dry with a dry towel. Wait until scab get rid of it, you can wash your face with facial cleanser.

### 2. How to take care

Can't drink before the scab falls off, and other can eat casually; If the area is large, and causes infection accidentally, you can buy Iodine to wipe, and then apply the professional medicines such as redlyticin ointment, Jingwanhong Scald Palace, etc., or follow the doctor's advice. In addition, let scab go off naturally, don't pull it with your hands!

### 3. About skin become red after using

There are generally red marks after crusting falling. Don't worry. This is common for skin growth. It is normal, that's the color of the new skin! Generally speaking, it will disappear by itself. If you want it to disappear faster, you can massage it. Take a small amount of vitamin E or aloe vera gel in the red seal. Use your fingers to rub it with your fingers.

### Method of operation

When used for the first time, you may have various problems and concerns. This is normal, mainly because it has not been operated. It is recommended to practice with a colorful thing. After practicing, start with your body arm, and finally operate on the face, step by step, so that you have confidence in yourself! If you do not understand, you can contact the customer service, and the customer service will guide you.

### Specifications

Model	LS-006
Size	209.4x41.4x37.0mm
Light source Type	laser
Wavelength	440~460nm
Rated power	< 8W
Adaptor	Output:AC100-240V, 50/60Hz Input:DC5V --- 2A
Work Tem.	Temperature: 5 ~ 40°C Relative humidity: 20 ~ 90%, no dew
Storage Tem.	Temperature: -20 ~ 55°C Relative humidity: 5 ~ 95%, no dew

## Simple fault exclude

Performance	Reason	Deal way
Cannot turn on	1. Press the "boot" button without force exertion; 2. Press and hold the "boot" not exceeded 2 seconds	When booting, please apply for a little force and hold the "Power button" for no less than 2s.
No laser shot when press shot button	Only press shot button will not have laser output	Correct operation: Press "Positioning button" first, don't remove, then press the "shot button"
Automatic shutdown	No operation in 5min.	Press "Boot" button to turn on



## Warranty Card

### Dear "MOOLWEEL" product user:

First of all, thank you very much for choosing "Moolweel" products, so that we have the opportunity to provide you with the best quality service. In order to make our services more satisfied, please read this instructions carefully and keep this warranty properly after buying Card.

### Warranty

During the warranty period, the protection of the Moolweel products (for free maintenance) is implemented from the date of purchase. The products during the warranty period are only maintained at the brand identification location; during warranty, please show this warranty card, purchase time and voucher. During the date of purchase, during the following period, it is repaired under normal use; the whole machine warranty (except commercial products) is one of the following situations. Even during the warranty period, but it gave charges to repair.

- Product damage caused by violating the instructions of the product use; Product damage caused by unprecedented expansion components (such as chargers, etc.);
- Product damage caused by exceed the allowable use environment (such as high temperature, corrosion, etc.);
- Product damage caused by non-normal external forces (such as falling, squeezing, bumps, impacts, etc.);
- Natural disasters or other force majeure Products caused by earthquakes);

- Products that counterfeit and illegal smuggling; cannot provide valid purchase vouchers, shipments, warranty cards, etc.;
- Other products or parts caused by other unexpected factors or useless users;

In order to provide better service, pls fill in the below info.:

Name: \_\_\_\_\_

Wechat: \_\_\_\_\_

Buy ID: \_\_\_\_\_

Contact: \_\_\_\_\_

Warranty Reason: \_\_\_\_\_

装订方式：骑马订装

尺寸：54×120mm

材质：封面157g铜版纸

内页105g铜版纸

不用过膜

封二封三空白，单色印刷。