
40K SLIMMING MACHINE

< Instruction Manual >

In order to use this product correctly, please read the manual carefully and keep it properly.

Part name

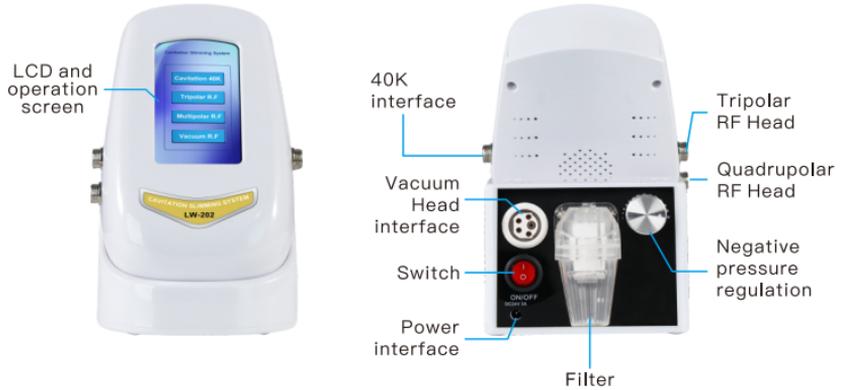Triolar
RF Head

Quadrupolar
RF Head

Vacuum Head


40K
Cavitation Head

Function of Different Head

Tripolar RF Head: Recommend level 1–2

Used for small areas, like eye area, face area

When used on the skin, it can firm and lift the skin, reduce wrinkles; improve eye circulation, shrink pores, and reduce puffiness; eliminate wrinkles, dark circles, and under-eye bags; enhances skin texture and eliminates sagging around the eyes; rejuvenate the skin and enhances eye contour.

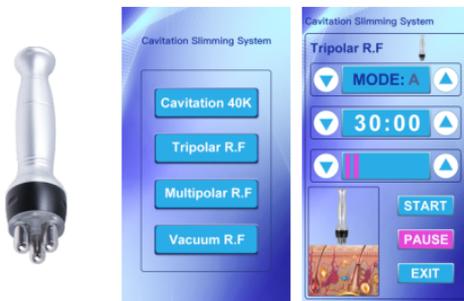When using, select Tripolar RF Head and enter the function panel.

Time: Default duration is set to 30 minutes, but it can be adjusted within a range of 0 to 60 minutes. If set to 0, the head will cease operation.

Energy: Higher the level, higher the energy.

There are 10 levels adjustable.

Press “START” to start the operation. Press the “PAUSE” to pause the operation. After use, press “EXIT” to back to

homepage.

MODE A: the red light flashes constantly

MODE B: the red light flashes at a frequency of 1 second intervals


MODE C: the red light flashes at a frequency of 2 second intervals

MODE D: the red light flashes at a frequency of 3 second intervals

Quadrupolar RF Head:Recommend level 2-3

Used for large area, like body area, arms, tummy, thigh

It helps firm the skin, enhance skin elasticity, eliminate facial wrinkles, tighten and lift the skin, reduce facial swelling and redness, relieve eye fatigue, remove dark circles and eye bags, and overall improve skin appearance.

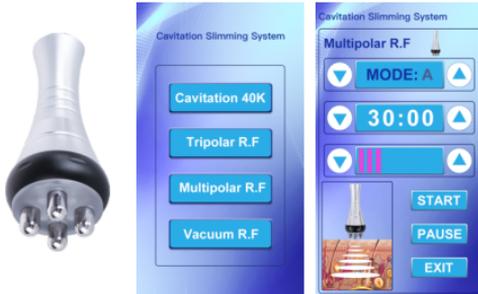When using, select Quadrupolar RF Head and enter the function panel.

Time: Default duration is set to 30 minutes, but it can be adjusted within a range of 0 to 60 minutes. If set to 0, the head will cease operation.

Energy: Higher the level, higher the energy.

There are 10 levels adjustable.

Press “START” to start the operation. Press the “PAUSE” to pause the operation. After use, press “EXIT” to back to homepage.

MODE A: the red light flashes constantly

MODE B: the red light flashes at a frequency of 1 second intervals

MODE C: the red light flashes at a frequency of 2 second intervals

MODE D: the red light flashes at a frequency of 3 second intervals

40k cavitation head

Used for body shaping, not face

It effectively dissolves fat, detoxifies lymph, tightens skin, enhances skin elasticity, and targets the fat layer directly. By stimulating sweat glands, liver intestine circulation, and lymph circulation, it efficiently eliminates excess fat and toxins from the body.

When using, select select 40k cavitation head and enter the function panel.

Time: Default duration is set to 30 minutes, but it can be adjusted within a range of 0 to 60 minutes. If set to 0, the head will cease operation.

Energy: Higher the level, higher the energy.

There are 10 levels adjustable.

Press “START” to start the operation. Press the “PAUSE” to pause the operation. After use, press “EXIT” to back to homepage.

MODE A: continuous vibration

MODE B: stopping for 1 second per vibration

MODE C: stopping for 2 seconds per vibration

MODE D: stopping for 3 seconds per vibration

Vacuum Head

Body maintenance

The diversified level adjustment ensures comfortable care for your body, while promoting collagen regeneration and elasticity to improve skin texture. It also addresses skin sagging caused by fat and aging, aids in waste elimination and metabolism promotion. The negative pressure massage function effectively reduces fat cell volume, allowing the multi-polar thermal melting radio frequency to penetrate deep into the skin tissue more efficiently. Simultaneously, the suction of negative pressure enhances blood circulation and metabolism, providing full relaxation for muscle tissue and relieving fatigue in both the skin and muscles.

Please make sure filter plugged in, otherwise no suction.

Twist the button to adjust vacuum intensity, default intensity is the maximum value. Twist the button counterclockwise to decrease the vacuum intensity. Twist the button clockwise to increase the vacuum intensity.

When using, select select vacuum head and enter the function panel.

Time: Default duration is set to 30 minutes, but it can be adjusted within a range of 0 to 60 minutes. If set to 0, the head will cease operation.

Energy: Adjust the energy level of RF.

Higher the level, higher the energy.

There are 10 levels adjustable.


(Please do not set suction time and release time same, it will not have suction)

Press “START” to start the operation. Press the “PAUSE” to pause the operation. After use, press “EXIT” to back to homepage.

Using Guidance


Thigh

- (1) From bottom to top, push to groin to dredge lymph;
- (2) Turning anticlockwise from bottom to top can help break down fat;
- (3) Decompose greasy by pushing back and forth;
- (4) Pull down from the knee to the groin to improve the curve.

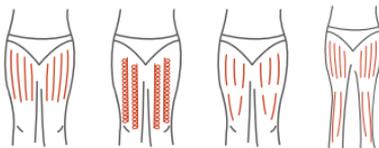
Hips

- (1) Pull to the waist in the direction of lymph;
- (2) Circle anticlockwise from bottom to top and move back and forth;
- (3) According to the muscle texture, it can slide and compact the buttock skin, making the buttock plump and firm, presenting the perfect curve.

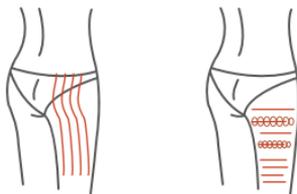
Back

- (1) Push back and forth 2–3 times;
- (2) Circle anticlockwise to stimulate the lymph nodes;
- (3) On the back, turn anticlockwise and move up and down to help break down fat;
- (4) Using lymphatic drainage technique to pull toxin to lymph nodes can help improve back curve.

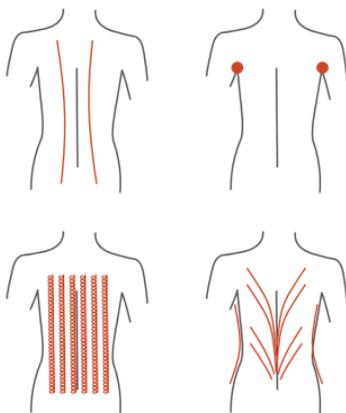Arm

- (1) From the inner elbow to the armpit, do lymphatic drainage;
- (2) To stimulate lymph nodes in a circle;
- (3) Break down excess fat by turning it anticlockwise or pushing it back and forth;
- (4) Tighten skin from elbow to armpit.

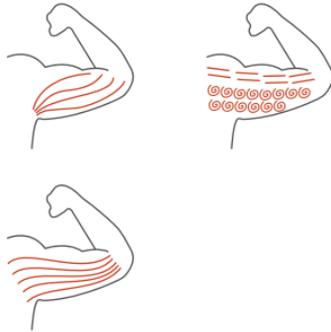
Abdomen

- (1) Around the navel from small to large, clockwise circle to help large intestine peristalsis;
- (2) From the small circle to the large circle, push it anticlockwise.

Face

- (1) Evenly apply the product on the face, use the face RF head, operate from the bottom to the top to lift the face, make the skin in a tense state, it is recommended to use for 5 minutes;
- (2) Apply the product evenly to the face, use the RF head on the face, operate in the direction of the arrow, it is recommended to use it for 5–8 minutes.

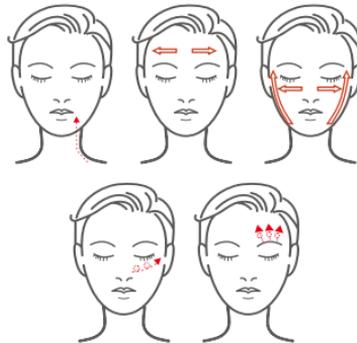Warm Tips for You

1. Do not use on your original and dry skin surface, it will burn your skin. Make sure you apply cream or oils to your skin then you can use.
2. When using, there must be an aqueous medium (essential oil/creams, etc.), otherwise it will damage to the machine Do not use on the same part for too long, and avoid body parts with bones.
3. When using, there must be an aqueous medium (essential oil/creams, etc.), otherwise it will damage to the machine. Do not connect the power until probe is well installed.
4. The machine is not allowed to replace probe when in the charging state.
5. Please take off the metal jewelry when using the machine.
6. When the machine is used for 1 hour, please stop for 10 minutes to let the machine rest.
7. Make sure you turn off the power after using, and clean the probe keep good maintenance for next time use.
8. Please keep the probe clean every time, use a dry cloth to wipe it, do not use water to clean.
9. After using, clean you skin and apply hydrating mask or creams, oils after 10 minutes.
10. When using 80KHead, some people may experience tinnitus. This is normal because this sound wave exceeds the range of human sound waves.

Safety Attentions

- Don't use on the skin area with surgery wound, scars wound.
- Don't use on the dermatology parts, eczema, and other sensitive skin.
- People who have hearing aids, heart pacemakers and heart machine.
- People who are minors, pregnant are prohibited to use.
- People who have heart disease, hypertension, coronary heart disease, diabetes, skin cancer, thyroid, malignancy, current sensitive group are prohibited to use.
- People who have infections of the skin, skin burns, anti-clotting disorders, skin bleeding, or skin allergies are prohibited to use.
- People with metal or silicone inside the body cannot use.

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